

Nanonis TCP Protocol

TCP Programming Interface

Table of Contents

General.....	12
Architecture	12
Data Types.....	12
Request message	14
Header.....	14
Body	14
Response message	14
Header.....	14
Body	14
Examples	15
Read limits of Step Channel 1 in 3D Sweeper	16
Change limits of Step Channel 1 in 3D Sweeper	18
Troubleshooting.....	21
Change Log.....	23
Functions.....	24
3D Sweeper	24
3DSwp.AcqChsSet	24
3DSwp.AcqChsGet.....	24
3DSwp.SaveOptionsSet.....	25
3DSwp.SaveOptionsGet	25
3DSwp.Start	26
3DSwp.Stop.....	26
3DSwp.Open	26
3DSwp.StatusGet	26
3DSwp.SwpChSignalSet.....	27
3DSwp.SwpChSignalGet.....	27
3DSwp.SwpChLimitsSet	27

3DSwp.SwpChLimitsGet	28
3DSwp.SwpChPropsSet	28
3DSwp.SwpChPropsGet	29
3DSwp.SwpChTimingSet	29
3DSwp.SwpChTimingGet	29
3DSwp.SwpChModeSet	30
3DSwp.SwpChModeGet	30
3DSwp.SwpChMLSSet	30
3DSwp.SwpChMLSGet	31
3DSwp.StpCh1SignalSet	31
3DSwp.StpCh1SignalGet	31
3DSwp.StpCh1LimitsSet	32
3DSwp.StpCh1LimitsGet	32
3DSwp.StpCh1PropsSet	32
3DSwp.StpCh1PropsGet	33
3DSwp.StpCh1TimingSet	33
3DSwp.StpCh1TimingGet	33
3DSwp.StpCh2SignalSet	34
3DSwp.StpCh2SignalGet	34
3DSwp.StpCh2LimitsSet	34
3DSwp.StpCh2LimitsGet	35
3DSwp.StpCh2PropsSet	35
3DSwp.StpCh2PropsGet	35
3DSwp.StpCh2TimingSet	36
3DSwp.StpCh2TimingGet	36
3DSwp.TimingRowLimitSet	36
3DSwp.TimingRowLimitGet	37
3DSwp.TimingRowMethodsSet	37
3DSwp.TimingRowMethodsGet	38
3DSwp.TimingRowValsSet	38
3DSwp.TimingRowValsGet	39
3DSwp.TimingEnable	39

3DSwp.TimingSend	39
3DSwp.FilePathsGet	40
1D Sweeper	41
1DSwp.AcqChsSet	41
1DSwp.AcqChsGet	41
1DSwp.SwpSignalSet	42
1DSwp.SwpSignalGet	42
1DSwp.LimitsSet	42
1DSwp.LimitsGet	43
1DSwp.PropsSet	43
1DSwp.PropsGet	43
1DSwp.Start	44
1DSwp.Stop	44
1DSwp.Open	44
1DSwp.ModeSet	45
1DSwp.ModeGet	45
1DSwp.SwpSignalListGet	46
High-Speed Sweeper	47
HSSwp.AcqChsSet	47
HSSwp.AcqChsGet	47
HSSwp.AutoReverseSet	47
HSSwp.AutoReverseGet	48
HSSwp.EndSettlSet	48
HSSwp.EndSettlGet	49
HSSwp.NumSweepsSet	49
HSSwp.NumSweepsGet	49
HSSwp.ResetSignalsSet	49
HSSwp.ResetSignalsGet	50
HSSwp.SaveBasenameSet	50
HSSwp.SaveBasenameGet	50
HSSwp.SaveDataSet	51
HSSwp.SaveDataGet	51

HSSwp.SaveOptionsSet	51
HSSwp.SaveOptionsGet	52
HSSwp.Start.....	52
HSSwp.Stop	52
HSSwp.StatusGet	53
HSSwp.SwpChSigListGet	53
HSSwp.SwpChSignalSet.....	53
HSSwp.SwpChSignalGet	54
HSSwp.SwpChLimitsSet.....	54
HSSwp.SwpChLimitsGet	54
HSSwp.SwpChNumPtsSet	55
HSSwp.SwpChNumPtsGet.....	55
HSSwp.SwpChTimingSet	55
HSSwp.SwpChTimingGet.....	56
HSSwp.SwpChBwdSwSet	56
HSSwp.SwpChBwdSwGet.....	56
HSSwp.SwpChBwdDelaySet	57
HSSwp.SwpChBwdDelayGet	57
HSSwp.StpChSigListGet	57
HSSwp.StpChSignalSet	58
HSSwp.StpChSignalGet.....	58
HSSwp.StpChLimitsSet	58
HSSwp.StpChLimitsGet	59
HSSwp.StpChNumPtsSet	59
HSSwp.StpChNumPtsGet	59
HSSwp.StpChTimingSet.....	59
HSSwp.StpChTimingGet	60
Signals	61
Signals.NamesGet	61
Signals.CalibrGet	61
Signals.RangeGet.....	61
Signals.ValGet	62

Signals.ValsGet	62
Signals.MeasNamesGet	63
Signals.AddRTGet	63
Signals.AddRTSet.....	64
User Inputs.....	65
UserIn.CalibrSet	65
User Outputs	66
UserOut.ModeSet	66
UserOut.ModeGet	66
UserOut.MonitorChSet	66
UserOut.MonitorChGet.....	67
UserOut.ValSet.....	67
UserOut.CalibrSet	67
UserOut.CalcSignalNameSet	68
UserOut.CalcSignalNameGet	68
UserOut.CalcSignalConfigSet	69
UserOut.CalcSignalConfigGet.....	70
UserOut.LimitsSet	71
UserOut.LimitsGet	71
UserOut.SlewRateSet.....	71
UserOut.SlewRateGet	72
Digital Lines	73
DigLines.PropsSet.....	73
DigLines.OutStatusSet.....	73
DigLines.TTLValGet.....	74
DigLines.Pulse	74
Data Logger	75
DataLog.Open	75
DataLog.Start	75
DataLog.Stop.....	75
DataLog.StatusGet	76
DataLog.ChsSet	76

DataLog.ChsGet.....	77
DataLog.PropsSet.....	77
DataLog.PropsGet.....	78
TCP Logger.....	79
TCPLog.Start.....	79
TCPLog.Stop.....	79
TCPLog.ChsSet.....	79
TCPLog.OversamplSet.....	80
TCPLog.StatusGet.....	80
Oscilloscope High Resolution.....	81
OsciHR.ChSet.....	81
OsciHR.ChGet.....	81
OsciHR.OversamplSet.....	81
OsciHR.OversamplGet.....	82
OsciHR.CalibrModeSet.....	82
OsciHR.CalibrModeGet.....	82
OsciHR.SamplesSet.....	83
OsciHR.SamplesGet.....	83
OsciHR.PreTrigSet.....	83
OsciHR.PreTrigGet.....	84
OsciHR.Run.....	84
OsciHR.OsciDataGet.....	84
OsciHR.TrigModeSet.....	85
OsciHR.TrigModeGet.....	85
OsciHR.TrigLevChSet.....	85
OsciHR.TrigLevChGet.....	86
OsciHR.TrigLevValSet.....	86
OsciHR.TrigLevValGet.....	86
OsciHR.TrigLevHystSet.....	86
OsciHR.TrigLevHystGet.....	87
OsciHR.TrigLevSlopeSet.....	87
OsciHR.TrigLevSlopeGet.....	87

OsciHR.TrigDigChSet	88
OsciHR.TrigDigChGet.....	88
OsciHR.TrigArmModeSet	88
OsciHR.TrigArmModeGet.....	89
OsciHR.TrigDigSlopeSet.....	89
OsciHR.TrigDigSlopeGet.....	89
OsciHR.TrigRearm	89
OsciHR.PSDShow	90
OsciHR.PSDWeightSet.....	90
OsciHR.PSDWeightGet	90
OsciHR.PSDWindowSet.....	91
OsciHR.PSDWindowGet	91
OsciHR.PSDAvgTypeSet	91
OsciHR.PSDAvgTypeGet.....	91
OsciHR.PSDAvgCountSet	92
OsciHR.PSDAvgCountGet.....	92
OsciHR.PSDAvgRestart.....	92
OsciHR.PSDDataGet	93
Script	94
Script.Open	94
Script.Load	94
Script.Save.....	94
Script.Deploy.....	95
Script.Undeploy.....	95
Script.Run	95
Script.Stop.....	95
Script.ChsGet.....	96
Script.ChsSet	96
Script.DataGet.....	97
Script.Autosave	97
Script.LUTOpen	98
Script.LUTLoad	98

Script.LUTSave.....	98
Script.LUTDeploy.....	99
Lock-In.....	100
LockIn.ModOnOffSet.....	100
LockIn.ModOnOffGet.....	100
LockIn.ModSignalSet.....	100
LockIn.ModSignalGet.....	101
LockIn.ModPhasRegSet.....	101
LockIn.ModPhasRegGet.....	102
LockIn.ModHarmonicSet.....	102
LockIn.ModHarmonicGet.....	103
LockIn.ModPhasSet.....	103
LockIn.ModPhasGet.....	103
LockIn.ModAmpSet.....	104
LockIn.ModAmpGet.....	104
LockIn.ModPhasFreqSet.....	104
LockIn.ModPhasFreqGet.....	105
LockIn.DemodSignalSet.....	105
LockIn.DemodSignalGet.....	105
LockIn.DemodHarmonicSet.....	106
LockIn.DemodHarmonicGet.....	106
LockIn.DemodHPFilterSet.....	107
LockIn.DemodHPFilterGet.....	107
LockIn.DemodLPFilterSet.....	108
LockIn.DemodLPFilterGet.....	108
LockIn.DemodPhasRegSet.....	109
LockIn.DemodPhasRegGet.....	109
LockIn.DemodPhasSet.....	110
LockIn.DemodPhasGet.....	110
LockIn.DemodSyncFilterSet.....	110
LockIn.DemodSyncFilterGet.....	111
LockIn.DemodRTSignalsSet.....	111

LockIn.DemodRTSignalsGet	112
Lock-In Frequency Sweep	113
LockInFreqSwp.Open	113
LockInFreqSwp.Start	113
LockInFreqSwp.SignalSet	114
LockInFreqSwp.SignalGet.....	114
LockInFreqSwp.LimitsSet	114
LockInFreqSwp.LimitsGet.....	114
LockInFreqSwp.PropsSet.....	115
LockInFreqSwp.PropsGet	115
Utilities	116
Util.SessionPathGet	116
Util.SessionPathSet	116
Util.SettingsLoad	116
Util.SettingsSave	117
Util.LayoutLoad	117
Util.LayoutSave	117
Util.Lock	118
Util.UnLock.....	118
Util.RTFreqSet	118
Util.RTFreqGet.....	118
Util.AcqPeriodSet	119
Util.AcqPeriodGet	119
Util.RTOversamplSet	119
Util.RTOversamplGet	119
Util.Quit.....	120
Util.VersionGet.....	120
File.....	121
File.datLoad.....	121
MCVA5 Preamplifier	122
MCVA5.UserInSet.....	122
MCVA5.UserInGet	122

MCVA5.GainSet.....	123
MCVA5.GainGet	123
MCVA5.InputModeSet.....	123
MCVA5.InputModeGet	124
MCVA5.CouplingSet	124
MCVA5.CouplingGet	124
MCVA5.SingleStateUpdate	125
MCVA5.ContStateUpdateSet	125
MCVA5.ContStateUpdateGet.....	125
MCVA5.SingleTempUpdate.....	126
MCVA5.ContTempUpdateSet	126
MCVA5.ContTempUpdateGet.....	126
Generic PI Controller.....	127
PICtrl.OnOffSet.....	127
PICtrl.OnOffGet	127
PICtrl.CtrlChSet.....	127
PICtrl.CtrlChGet.....	128
PICtrl.InputChSet.....	128
PICtrl.InputChGet	129
PICtrl.PropsSet	129
PICtrl.PropsGet.....	129
PICtrl.CtrlChPropsSet	130
PICtrl.CtrlChPropsGet.....	130

TCP Protocol

General

Architecture

The Nanonis software works as a TCP Server, and a remote application works as a TCP Client. The TCP Server listens at the ports specified in the Options window (under the System menu). The Client can open one, two, three or four different connections to the Server at these ports.

Each individual connection handles commands serially (one command after another) guaranteeing synchronized execution. These connections can be found and configured in the Options window (under the System menu).

Every message sent from the client to the server (request message), and viceversa (response message) when *Send response back* is set to True in the request message (see Request message>Header section), consists of header and body.

All numeric values are sent in binary form (e.g. a 32 bit integer is encoded in 4 bytes). The storage method of binary encoded numbers is big-endian, that is, the most significant byte is stored at the lowest address.

Data Types

There are thirteen data types which appear in the header and body of both request message and response message:

string

This is an array of characters, where every character has a size of one byte. In the header, the strings have a fixed size (*Command name* is 32 bytes), and in the body, the strings have a variable size which is always prepended as an integer 32.

int

32 bit signed integer. Its range is -2147483648 to 2147483647.

unsigned int16

16 bit unsigned integer. Its range is 0 to 65535.

unsigned int32

32 bit unsigned integer. Its range is 0 to 4294967295.

float32

32 bit (single precision) floating point number.

float64

64 bit (double precision) floating point number.

1D array string

This is a unidimensional array of strings. The array size (number of elements) and its size in bytes are sent as independent arguments right before the array. Each element of the array is obviously a string, preceded by its size. See Functions for more details.

1D array int

This is a unidimensional array of 32 bit signed integers. The array size (number of elements) is usually sent as an independent argument right before the array. See Functions for more details.

1D array unsigned int8

This is a unidimensional array of 8 bit unsigned integers. The array size (number of elements) is usually sent as an independent argument right before the array. See Functions for more details.

1D array unsigned int32

This is a unidimensional array of 32 bit unsigned integers. The array size (number of elements) is usually sent as an independent argument right before the array. See Functions for more details.

1D array float32

This is a unidimensional array of 32 bit floating point numbers. The array size (number of elements) is usually sent as an independent argument right before the array. See Functions for more details.

1D array float64

This is a unidimensional array of 64 bit floating point numbers. The array size (number of elements) is usually sent as an independent argument right before the array. See Functions for more details.

2D array float32

This is a bi-dimensional array of 32 bit floating point numbers. The array size (number of rows and columns) is usually sent as two independent arguments right before the array. See Functions for more details.

2D array string

This is a bi-dimensional array of strings. The array size (number of rows and columns) is usually sent as two independent arguments right before the array. Each element of the array is obviously a string, preceded by its size. See Functions for more details

Request message

This is the message sent from the client to the server.

Header

Header size is fixed to 40 bytes (last 2 bytes are currently not used and they should be set to zero) and contains the following elements:

- **Command name** (string) (32 bytes) is the name of the executed command. It matches one of the function names described in the Functions section of this document (i.e. *BiasSpectr.Open*). Maximum number of characters is 32.
- **Body size** (int) (4 bytes) is the size of the message body in bytes.
- **Send response back** (unsigned int16) (2 bytes) defines if the server sends a message back to the client (=1) or not (=0). All functions can return at least the error information returned after executing the specified function.

Body

The body size is variable and contains the argument values (if any) sent to the server. Each function has its own arguments described in detail in the Functions section.

Response message

This is only sent from the server to the client if the flag to send the response back (in the request message to the server) is true.

Be aware that without the response message the order of execution between different TCP connections cannot be guaranteed. On the other hand, the order of execution of commands sent through the same TCP connection is guaranteed as the commands are serialized.

Header

Header size is fixed to 40 bytes (last 4 bytes are currently not used and they should be set to zero) and contains the following elements:

- **Command name** (string) (32 bytes) is the name of the executed command. It matches one of the function names described in the Functions section of this document (i.e. *BiasSpectr.Open*). Maximum number of characters is 32.
- **Body size** (int) (4 bytes) is the size of the message body in bytes.

Body

The body size is variable and contains the argument values returned by the function (sent from the server to the client). Each function has its own return arguments described in detail in the Functions section.

After the return arguments values, the body includes the error information containing the following elements:

- **Error status** (unsigned int32) (4 bytes) returns 1=True if there is an error when executing the function.
- **Error description size** (int) (4 bytes) returns the size of the error description which follows.
- **Error description** (string) (variable size) returns the description of the error.

Examples

The following examples show the encoded strings sent through TCP/IP in the request message to execute the desired functions, i.e. from the client (remote application) to the server (Nanonis software).

It is also explained the encoded strings in the response message when the *Send response back* flag in the request message is set to True, i.e. received by the client (remote application) from the server (Nanonis software).

The commands are displayed in hexadecimal (one hexadecimal number pair corresponds to one byte):

Read limits of Step Channel 1 in 3D Sweeper

This example reads the limits of Step Channel 1 in the 3D Sweeper module. It uses the function *3DSwp.StpCh1LimitsGet*.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops
- **Error** described in the Response message>Body section

REQUEST MESSAGE:

The Header is fixed always to 40 bytes (last 2 bytes not used), containing the following elements:

- **Command name (string) (32 bytes)** is the hexadecimal representation of the command name *3DSwp.StpCh1LimitsGet* padded with zeros to length 32:
3364 7377 702E 7374 7063 6831 6C69 6D69 7473 6765 7400 0000 0000 0000 0000
- **Body size (int) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to the integer size of the message body in bytes (=0 because there are no arguments) padded with zeros to length 4:
0000 0000
- **Send response back (unsigned int16) (2 bytes)** defines if the server sends a message back to the client. In this example we set it to True (=1):
0001
- **Not used (2 bytes):**
0000

The Body is of variable size and contains the argument values (if any) sent to the server. In this example this function has no arguments in the Request message.

Request message					
	Header				Body
	Command name	Body size in bytes	Send response back	Not Used	-
Size (Bytes)	Fixed (32)	Fixed (4)	Fixed (2)	Fixed (2)	-
Readable value representation	3DSwp.StpCh1LimitsGet	0	True	-	-
Hex representation of string to be sent over TCP	3364 7377 702E 7374 7063 6831 6C69 6D69 7473 6765 7400 0000 0000 0000 0000	0000 0000	0001	0000	-

RESPONSE MESSAGE:

The Header is fixed always to 40 bytes (last 4 bytes not used), containing the following elements:

- **Command name (string) (32 bytes)** is the hexadecimal representation of the command name *3DSwp.StpChlLimitsGet* padded with zeros to length 32 (like in the request message):
3364 7377 702E 7374 7063 6831 6C69 6D69 7473 6765 7400 0000 0000 0000 0000
- **Body size (int) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to the integer size of the message body in bytes (=16 bytes of the arguments *Start*, *Stop*, *Error status*, and *Error size*) padded with zeros to length 4:
0000 0010
- **Not used (4 bytes):**
0000 0000

The Body is of variable size and contains the argument values (if any) sent from the server to the client. In this example this function has the following arguments:

- **Start (float32) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to -1mm:
BA83 126F
- **Stop (float32) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to 2mm:
3B03 126F
- **Error status (unsigned int32) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to 0 (no error):
0000 0000
- **Error size (int) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to the integer size of the error description in bytes (=0 in this case because there is no error):
0000 0000
- **Error description (string) (variable)** is the hexadecimal representation of the big-endian encoded string corresponding to the integer size of the error description in bytes (=0 in this case because there is no error):
0000 0000

Response message								
	Header			Body				
	Command name	Body size in bytes	Not Used	Start	Stop	Error status	Error size in bytes	Error description
Size (Bytes)	Fixed (32)	Fixed (4)	Fixed (4)	4	4	4	4	0
Readable value representation	3DSwp.StpChlLimitsGet	16	-	-1m	2m	False	0	-
Hex representation of string to be sent over TCP	3364 7377 702E 7374 7063 6831 6C69 6D69 7473 6765 7400 0000 0000 0000 0000 0000	0000 0010	0000 0000	BA83 126F	3B03 126F	0000 0000	0000 0000	-

Change limits of Step Channel 1 in 3D Sweeper

This example changes the limits of Step Channel 1 in the 3D Sweeper module. It uses the function *3DSwp.StpCh1LimitsSet*.

Arguments:

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

REQUEST MESSAGE:

The Header is fixed always to 40 bytes (last 2 bytes not used), containing the following elements:

- **Command name (string) (32 bytes)** is the hexadecimal representation of the command name *3DSwp.StpCh1LimitsSet* padded with zeros to length 32:
3364 7377 702E 7374 7063 6831 6C69 6D69 7473 7365 7400 0000 0000 0000 0000
- **Body size (int) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to the integer size of the message body in bytes (=8 bytes of the arguments *Start* and *Stop*) padded with zeros to length 4:
0000 0008
- **Send response back (unsigned int16) (2 bytes)** defines if the server sends a message back to the client. In this example we set it to True (=1):
0001
- **Not used (2 bytes):**
0000

The Body is of variable size and contains the argument values (if any) sent to the server. In this example this function has the following arguments:

- **Start (float32) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to -3mm:
BB44 9BA6
- **Stop (float32) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to 4mm:
3B83 126F

Request message						
	Header				Body	
	Command name	Body size in bytes	Send response back	Not Used	Start	Stop
Size (Bytes)	Fixed (32)	Fixed (4)	Fixed (2)	Fixed (2)	-3m	4m
Readable value representation	3DSwp.StpCh1LimitsSet	8	True	-	10n	15n
Hex representation of string to be sent over TCP	3364 7377 702E 7374 7063 6831 6C69 6D69 7473 7365 7400 0000 0000 0000 0000 0000	0000 0008	0001	0000	BB44 9BA6	3B83 126F

RESPONSE MESSAGE:

The Header is fixed always to 40 bytes (last 4 bytes not used), containing the following elements:

- **Command name (string) (32 bytes)** is the hexadecimal representation of the command name *3DSwp.StpCh1LimitsSet* padded with zeros to length 32 (like in the request message):
3364 7377 702E 7374 7063 6831 6C69 6D69 7473 7365 7400 0000 0000 0000 0000
- **Body size (int) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to the integer size of the message body in bytes (=8 bytes of the arguments *Error status*, and *Error size*) padded with zeros to length 4:
0000 0008
- **Not used (4 bytes):**
0000 0000

The Body is of variable size and contains the argument values (if any) sent from the server to the client. In this example this function has the following arguments:


- **Error status (unsigned int32) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to 0 (no error):
0000 0000
- **Error size (int) (4 bytes)** is the hexadecimal representation of the big-endian encoded string corresponding to the integer size of the error description in bytes (=0 in this case because there is no error):
0000 0000
- **Error description (string) (variable)** is the hexadecimal representation of the big-endian encoded string corresponding to the integer size of the error description in bytes (=0 in this case because there is no error):

Response message						
	Header			Body		
	Command name	Body size in bytes	Not Used	Error status	Error size in bytes	Error description
Size (Bytes)	Fixed (32)	Fixed (4)	Fixed (4)	4	4	0
Readable value representation	3DSwp.StpCh1LimitsSet	8	-	False	0	-
Hex representation of string to be sent over TCP	3364 7377 702E 7374 7063 6831 6C69 6D69 7473 7365 7400 0000 0000 0000 0000	0000 0008	0000 0000	0000 0000	0000 0000	-

Troubleshooting

- 1) Make sure that the TCP connections between the TCP client and the Nanonis software (TCP server) are successfully established by using the corresponding remote ports specified in the Nanonis. The default ports can be changed.


This is available in the TCP Programming Interface section of the Main Options under the System menu. The displayed message shows if the connection has been established:

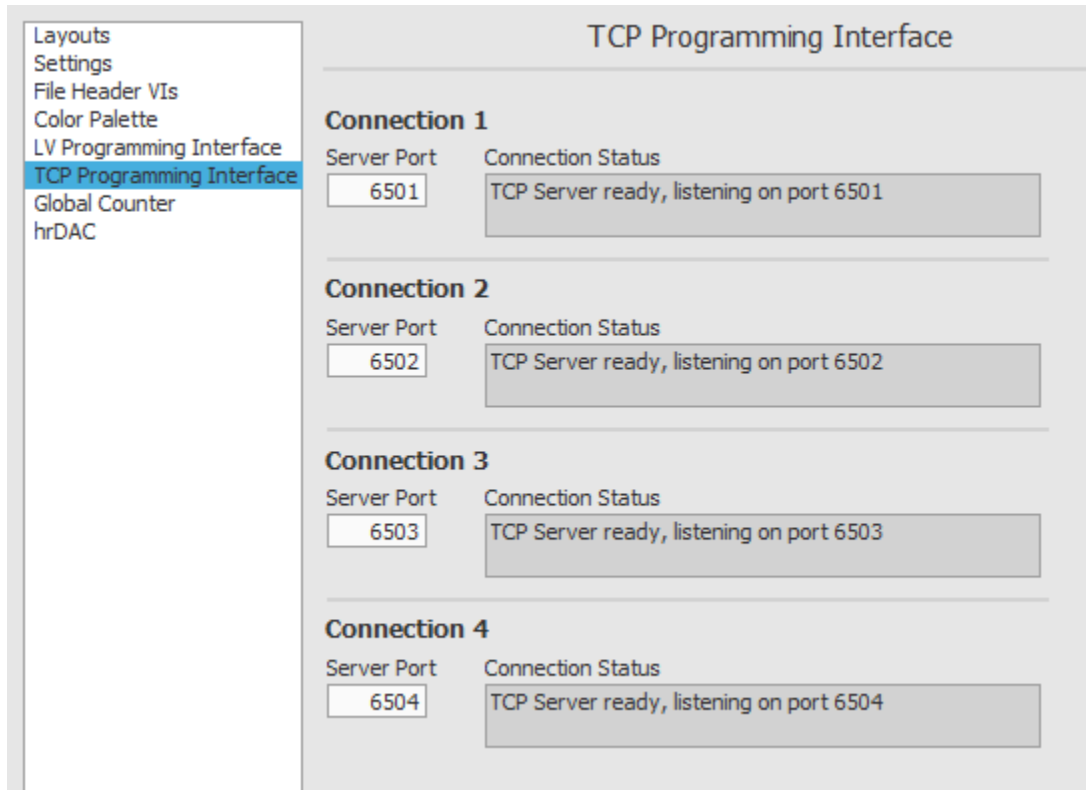

The screenshot displays the 'TCP Programming Interface' window. On the left is a sidebar menu with the following items: Layouts, Settings, File Header VIs, Color Palette, LV Programming Interface, TCP Programming Interface (highlighted in blue), Global Counter, and hrDAC. The main area of the window is titled 'TCP Programming Interface' and contains four sections, each for a different connection:

- Connection 1:** Server Port is 6501. Connection Status is 'TCP Server ready, listening on port 6501'.
- Connection 2:** Server Port is 6502. Connection Status is 'TCP Server ready, listening on port 6502'.
- Connection 3:** Server Port is 6503. Connection Status is 'TCP Server ready, listening on port 6503'.
- Connection 4:** Server Port is 6504. Connection Status is 'TCP Server ready, listening on port 6504'.

Before a TCP Connection is established

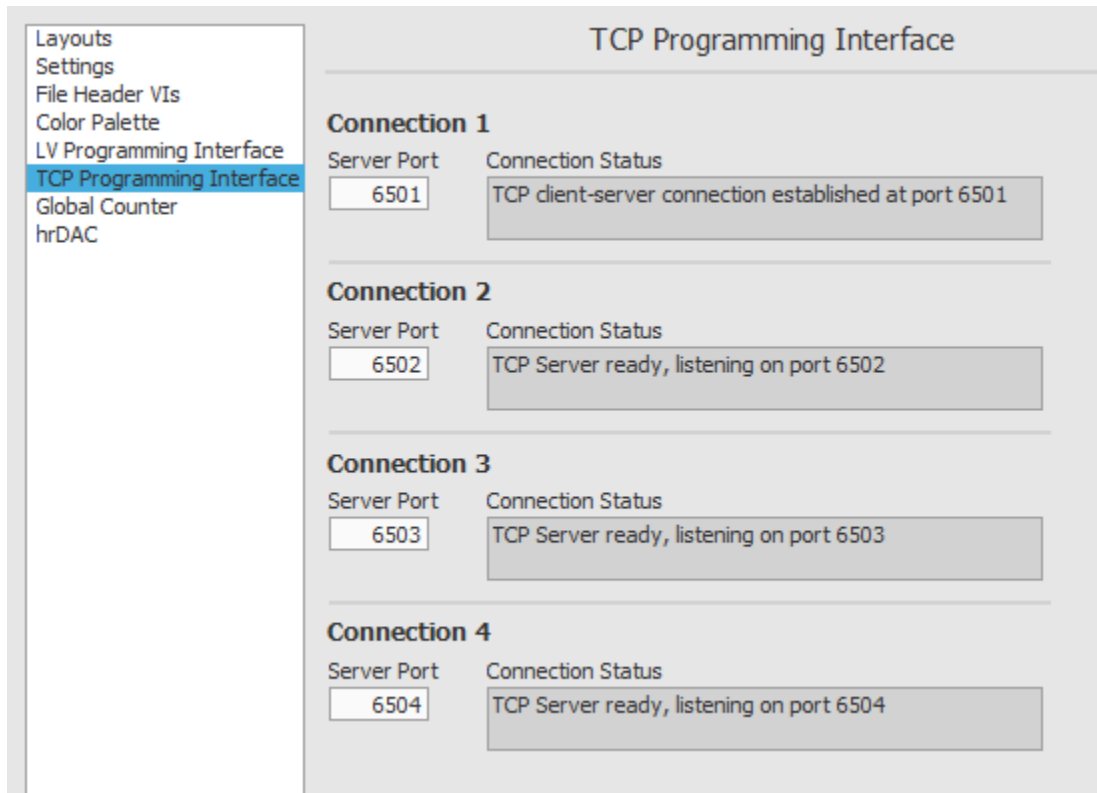

After a TCP Connection is established

- 2) When sending the TCP commands, the strings (like the command names) should NOT be sent using their hexadecimal representation.

In the examples we use the hexadecimal representation of a string as a way to explain the functions because a string might contain non-printable characters.

- 3) In the Request message, set the “Send response back” always to true to get a Response message from the Nanonis software.

If the formatting of the Request message is correct, you will get at least the error information (if any) when executing the function in the Nanonis software.

Change Log

R5.1.14766

Some functions from the Oscilloscope High Resolution have been adapted to include channel selection to support the 4-Channel version of this graph.

R5.1.14219

Added support in the Digital Lines functions for the expanded DIO ports available in selected systems through a dedicated DIO card.

Added functions to control the Generic PI Controllers.

R5.1.14161

Added functions to control the MCVA5 preamplifiers.

Added the function to switch the mode in the 1D Sweeper.

Added functions to control LUTs in the script module.

Added functions to open the script module and the LUT editor.

Changed the script load-save functions to add the script index.

Changed the data autosave function to add the saved path and basename.

Functions

3D Sweeper

3DSwp.AcqChsSet

Sets the list of recorded channels of the 3D Sweeper.

Arguments:

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of recorded channels. The indexes are comprised between 0 and 127, and it corresponds to the full list of signals available in the system.
To get the signal name and its corresponding index in the list of the 128 available signals in the Nanonis Controller, use the *Signal.NamesGet* function, or check the RT Idx value in the Signals Manager module.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.AcqChsGet

Returns the list of recorded channels of the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of the recorded channels. The indexes are comprised between 0 and 127, and it corresponds to the full list of signals available in the system
- **Channels names size** (int) is the size in bytes of the channels names array
- **Channels names number** (int) is the number of elements of the channels names array
- **Channels names** (1D array string) returns an array of channel names strings, where each string comes prepended by its size in bytes
- **Error** described in the Response message>Body section

3DSwp.SaveOptionsSet

Sets the saving options of the 3D Sweeper.

Arguments:

- **Series name size** (int) is the size (number of characters) of the series name string
- **Series name** (string) is the base name used for the saved sweeps. If empty string, there is no change
- **Create Date&Time Folder** (int) defines if this feature is active, where 0=no change, 1=On, 2=Off. If On, it creates a subfolder within the Session folder whose name is a combination of the basename and current date&time of the sweep, every time a sweep finishes.
- **Comment size** (int) is the size (number of characters) of the comment string
- **Comment** (string) is the comment saved in the header of the files. If empty string, there is no change
- **Modules names size** (int) is the size in bytes of the modules array. These are the modules whose parameters are saved in the header of the files
- **Modules names number** (int) is the number of elements of the modules names array
- **Modules names** (1D array string) is an array of modules names strings, where each string comes prepended by its size in bytes

Return arguments (if Send response back flag is set to True when sending request message to the server):

- **Error** described in the Response message>Body section

3DSwp.SaveOptionsGet

Returns the saving options of the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message to the server):

- **Series name size** (int) is the size of the series name string
- **Series name** (string) is the base name used for the saved sweeps
- **Create Date&Time Folder** (unsigned int32) returns if this feature is active, where 0=Off, 1=On
- **Fixed parameters size** (int) is the size in bytes of the Fixed parameters string array
- **Number of fixed parameters** (int) is the number of elements of the Fixed parameters string array
- **Fixed parameters** (1D array string) returns the fixed parameters of the sweep. The size of each string item comes right before it as integer 32.
- **Comment size** (int) is the size (number of characters) of the comment string
- **Comment** (string) is the comment saved in the header of the files
- **Modules parameters size** (int) is the size in bytes of the modules parameters array. These are the modules parameters saved in the header of the files
- **Modules parameters number** (int) is the number of elements of the modules parameters array
- **Modules parameters** (1D array string) is an array of modules names strings, where each string comes prepended by its size in bytes.
- **Error** described in the Response message>Body section

3DSwp.Start

Starts a sweep in the 3D Sweeper module.

When Send response back is set to True, it returns immediately afterwards.

Arguments:

- **Wait until finished** (unsigned int32) determines if the function returns immediately after starting a sweep (=0), or if the function wait until the sweep finishes (=1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.Stop

Stops the sweep in the 3D Sweeper module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.Open

Opens the 3D Sweeper module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.StatusGet

Returns the status of the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Status** (unsigned int32) is status of the 3D Sweep, where 0=Stopped, 1=Running, 2=Paused
- **Error** described in the Response message>Body section

3DSwp.SwpChSignalSet

Sets the Sweep Channel signal in the 3D Sweeper.

Arguments:

- **Sweep channel index** (int) is the index of the Sweep Channel, where -1 sets the *Unused* option

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.SwpChSignalGet

Returns the selected Sweep Channel signal in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Sweep channel index** (int) is the index of the Sweep Channel, where -1 is the *Unused* option
- **Channels names size** (int) is the size in bytes of the Channels names string array
- **Number of channels** (int) defines the number of elements of the Channels names string array
- **Channels names** (1D array string) returns the list of channels names. The size of each string item comes right before it as integer 32.
- **Error** described in the Response message>Body section

3DSwp.SwpChLimitsSet

Sets the limits of the Sweep Channel in the 3D Sweeper.

Arguments:

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.SwpChLimitsGet

Returns the limits of the Sweep Channel in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops
- **Error** described in the Response message>Body section

3DSwp.SwpChPropsSet

Sets the configuration of the Sweep Channel parameters in the 3D Sweeper.

Arguments:

- **Number of points** (int) sets the number of points of the sweep. 0 points means no change
- **Number of sweeps** (int) sets the total number of sweeps. 0 sweeps means no change
- **Backward sweep** (int) defines if the backward sweep is active, where 0=no change, 1=On, 2=Off
- **End of sweep action** (int) defines the behavior of the signal at the end of the sweep, where 0=no change, 1=no action, 2=reset signal to the original value, 3=go to arbitrary value
- **End of sweep arbitrary value** (float32) sets the arbitrary value to go at the end of the sweep if Go to arbitrary value is configured
- **Save all** (int) defines if all the configured sweeps are saved or only the averaged sweep, where 0=no change, 1=On, 2=Off

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.SwpChPropsGet

Returns the configuration of the Sweep Channel parameters in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of points** (int) returns the number of points of the sweep
- **Number of sweeps** (int) returns the total number of sweeps
- **Backward sweep** (unsigned int32) returns if the backward sweep is active, where 0=Off, 1=On
- **End of sweep action** (unsigned int32) returns the behavior of the signal at the end of the sweep, where 0=no action, 1=reset signal to the original value, 2=go to arbitrary value
- **End of sweep arbitrary value** (float32) returns the arbitrary value to go at the end of the sweep if Go to arbitrary value is configured
- **Save all** (unsigned int32) returns if all the configured sweeps are saved or only the averaged sweep, where 0=Off, 1=On
- **Error** described in the Response message>Body section

3DSwp.SwpChTimingSet

Sets the timing parameters of the Sweep Channel in the 3D Sweeper.

Arguments:

- **Initial settling time (s)** (float32)
- **Settling time (s)** (float32)
- **Integration time (s)** (float32)
- **End settling time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.SwpChTimingGet

Returns the timing parameters of the Sweep Channel in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Initial settling time (s)** (float32)
- **Settling time (s)** (float32)
- **Integration time (s)** (float32)
- **End settling time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)
- **Error** described in the Response message>Body section

3DSwp.SwpChModeSet

Sets the segments mode of the Sweep Channel signal in the 3D Sweeper.

Arguments:

- **Segments mode** (int) is the number of characters of the segments mode string.
If the segments mode is *Linear*, this value is 6. If the segments mode is *MLS*, this value is 3
- **Segments mode** (string) is *Linear* in Linear mode or *MLS* in MultiSegment mode

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.SwpChModeGet

Returns the segments mode of the Sweep Channel signal in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Segments mode** (int) is the number of characters of the segments mode string.
If the segments mode is *Linear*, this value is 6. If the segments mode is *MLS*, this value is 3
- **Segments mode** (string) is *Linear* in Linear mode or *MLS* in MultiSegment mode
- **Error** described in the Response message>Body section

3DSwp.SwpChMLSSet

Sets the MultiSegment values of the Sweep Channel in the 3D Sweeper.

Arguments:

- **Number of segments** (int) is the total number of segments. It defines the size of the following arrays
- **Segment Start values** (1D array float32) are the start values of the segments of the Sweep Channel
- **Segment Stop values** (1D array float32) are the stop values of the segments of the Sweep Channel
- **Segment Settling times** (1D array float32) are the settling times of the segments in seconds
- **Segment Integration times** (1D array float32) are the integration times of the segments in seconds
- **Segment Number of steps** (1D array int) are the number of steps of each segment
- **Last segment? array** (1D array unsigned int32) defines if the segments are the last one (1) or not (0)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.SwpChMLSGet

Returns the MultiSegment values of the Sweep Channel in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of segments** (int) is the total number of segments. It defines the size of the following arrays
- **Segment Start values** (1D array float32) are the start values of the segments of the Sweep Channel
- **Segment Stop values** (1D array float32) are the stop values of the segments of the Sweep Channel
- **Segment Settling times** (1D array float32) are the settling times of the segments in seconds
- **Segment Integration times** (1D array float32) are the integration times of the segments in seconds
- **Segment Number of steps** (1D array int) are the number of steps of each segment
- **Last segment? array** (1D array unsigned int32) defines if the segments are the last one (1) or not (0)
- **Error** described in the Response message>Body section

3DSwp.StpCh1SignalSet

Sets the Step Channel 1 signal in the 3D Sweeper.

Arguments:

- **Step channel 1 index** (int) is the index of the Step Channel 1, where -1 sets the *Unused* option

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.StpCh1SignalGet

Returns the selected Step Channel 1 signal in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Step channel 1 index** (int) is index of the Step Channel 1, where -1 is the *Unused* option
- **Channels names size** (int) is the size in bytes of the Channels names string array
- **Number of channels** (int) is the number of elements of the Channels names string array
- **Channels names** (1D array string) returns the list of channels names. The size of each string item comes right before it as integer 32.
- **Error** described in the Response message>Body section

3DSwp.StpCh1LimitsSet

Sets the limits of the Step Channel 1 in the 3D Sweeper.

Arguments:

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.StpCh1LimitsGet

Returns the limits of the Step Channel 1 in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops
- **Error** described in the Response message>Body section

3DSwp.StpCh1PropsSet

Sets the configuration of the Step Channel 1 parameters in the 3D Sweeper.

Arguments:

- **Number of points** (int) sets the number of points of the sweep. 0 points means no change
- **Backward sweep** (int) defines if the backward sweep is active, where 0=no change, 1=On, 2=Off
- **End of sweep action** (int) defines the behavior of the signal at the end of the sweep, where 0=no change, 1=no action, 2=reset signal to the original value, 3=go to arbitrary value
- **End of sweep arbitrary value** (float32) sets the arbitrary value to go at the end of the sweep if Go to arbitrary value is configured

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.StpCh1PropsGet

Returns the configuration of the Step Channel 1 parameters in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of points** (int) returns the number of points of the sweep
- **Backward sweep** (unsigned int32) returns if the backward sweep is active, where 0=Off, 1=On
- **End of sweep action** (unsigned int32) returns the behavior of the signal at the end of the sweep, where 0=no action, 1=reset signal to the original value, 2=go to arbitrary value
- **End of sweep arbitrary value** (float32) returns the arbitrary value to go at the end of the sweep if Go to arbitrary value is configured
- **Error** described in the Response message>Body section

3DSwp.StpCh1TimingSet

Sets the timing parameters of the Step Channel 1 in the 3D Sweeper.

Arguments:

- **Initial settling time (s)** (float32)
- **End settling time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.StpCh1TimingGet

Returns the timing parameters of the Step Channel 1 in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Initial settling time (s)** (float32)
- **End settling time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)
- **Error** described in the Response message>Body section

3DSwp.StpCh2SignalSet

Sets the Step Channel 2 signal in the 3D Sweeper.

Arguments:

- **Step Channel 2 name size** (int) is the number of characters of the Step Channel 2 name string
- **Step Channel 2 name** (string) is the name of the signal selected for the Step Channel 2

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.StpCh2SignalGet

Returns the selected Step Channel 2 signal in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Step Channel 2 name size** (int) is the number of characters of the Step Channel 2 name string
- **Step Channel 2 name** (string) is the name of the signal selected for the Step Channel 2
- **Channels names size** (int) is the size in bytes of the Channels names string array
- **Number of channels** (int) is the number of elements of the Channels names string array
- **Channels names** (1D array string) returns the list of channels names. The size of each string item comes right before it as integer 32.
- **Error** described in the Response message>Body section

3DSwp.StpCh2LimitsSet

Sets the limits of the Step Channel 2 in the 3D Sweeper.

Arguments:

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.StpCh2LimitsGet

Returns the limits of the Step Channel 2 in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops
- **Error** described in the Response message>Body section

3DSwp.StpCh2PropsSet

Sets the configuration of the Step Channel 2 parameters in the 3D Sweeper.

Arguments:

- **Number of points** (int) sets the number of points of the sweep. 0 points means no change
- **Backward sweep** (int) defines if the backward sweep is active, where 0=no change, 1=On, 2=Off
- **End of sweep action** (int) defines the behavior of the signal at the end of the sweep, where 0=no change, 1=no action, 2=reset signal to the original value, 3=go to arbitrary value
- **End of sweep arbitrary value** (float32) sets the arbitrary value to go at the end of the sweep if Go to arbitrary value is configured

Return arguments (if Send response back flag is set to True when sending request message):

Error described in the Response message>Body section

3DSwp.StpCh2PropsGet

Returns the configuration of the Step Channel 2 parameters in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of points** (int) returns the number of points of the sweep
- **Backward sweep** (unsigned int32) returns if the backward sweep is active, where 0=Off, 1=On
- **End of sweep action** (unsigned int32) returns the behavior of the signal at the end of the sweep, where 0=no action, 1=reset signal to the original value, 2=go to arbitrary value
- **End of sweep arbitrary value** (float32) returns the arbitrary value to go at the end of the sweep if Go to arbitrary value is configured
- **Error** described in the Response message>Body section

3DSwp.StpCh2TimingSet

Sets the timing parameters of the Step Channel 2 in the 3D Sweeper.

Arguments:

- **Initial settling time (s)** (float32)
- **End settling time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.StpCh2TimingGet

Returns the timing parameters of the Step Channel 2 in the 3D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Initial settling time (s)** (float32)
- **End settling time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)
- **Error** described in the Response message>Body section

3DSwp.TimingRowLimitSet

Sets the maximum time (seconds) and channel of the selected row in the Advanced Timing section of the 3D Sweeper.

Arguments:

- **Row index** (int) starting from 0 index
- **Maximum time (seconds)** (float32) defines the ultimate stop condition which is required since certain signal types can result in a target SNR or StdDev never being reached (infinite integration). Setting it to the minimum essentially switches off that set (limits to a single RT Cycle). NaN means no change
- **Channel index** (int) defines the channel to which the advanced configuration of the selected row is applied. -1 means no change

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.TimingRowLimitGet

Returns the maximum time (seconds) and channel of the selected row in the Advanced Timing section of the 3D Sweeper.

Arguments:

- **Row index** (int) starting from 0 index

Return arguments (if Send response back flag is set to True when sending request message):

- **Maximum time (seconds)** (float32) defines the ultimate stop condition which is required since certain signal types can result in a target SNR or StdDev never being reached (infinite integration).
- Setting it to the minimum essentially switches off that set (limits to a single RT Cycle).
- **Channel index** (int) defines the channel to which the advanced configuration of the selected row is applied
- **Channels names size** (int) is the size in bytes of the Channels names string array
- **Number of channels** (int) is the number of elements of the Channels names string array
- **Channels names** (1D array string) returns the list of channels names. The size of each string item comes right before it as integer 32.
- **Error** described in the Response message>Body section

3DSwp.TimingRowMethodsSet

Sets the methods of the selected row in the Advanced Timing section of the 3D Sweeper.

The possible values are -1=no change, 0=None, 1=Time, 2=Standard Deviation, 3=Signal to Noise Ratio.

Arguments:

- **Row index** (int) starting from 0 index
- **Method lower** (int) defines the method in the lower range of the selected row
- **Method middle** (int) defines the method in the middle range of the selected row
- **Method upper** (int) defines the method in the upper range of the selected row
- **Method alternative** (int) defines the method in the alternative range of the selected row

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.TimingRowMethodsGet

Returns the methods of the selected row in the Advanced Timing section of the 3D Sweeper.

The possible values are 0=None, 1=Time, 2=Standard Deviation, 3=Signal to Noise Ratio.

Arguments:

- **Row index** (int) starting from 0 index

Return arguments (if Send response back flag is set to True when sending request message):

- **Method lower** (int) gets the method in the lower range of the selected row
- **Method middle** (int) gets the method in the middle range of the selected row
- **Method upper** (int) gets the method in the upper range of the selected row
- **Method alternative** (int) gets the method in the alternative range of the selected row
- **Error** described in the Response message>Body section

3DSwp.TimingRowValsSet

Sets the ranges of the selected row in the Advanced Timing section of the 3D Sweeper.

Arguments:

- **Row index** (int) starting from 0 index
- **Middle range: from** (float64) is the upper limit of the lower range of the selected row
- **Lower range: value** (float64) is the value in the lower range of the selected row
- **Middle range: value** (float64) is the value in the middle range of the selected row
- **Middle range: to** (float64) is the lower limit of the upper range of the selected row
- **Upper range: value** (float64) is the value of the upper range of the selected row
- **Alternative range: value** (float64) is the value of the alternative range of the selected row

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.TimingRowValsGet

Returns the ranges of the selected row in the Advanced Timing section of the 3D Sweeper.

Arguments:

- **Row index** (int) starting from 0 index

Return arguments (if Send response back flag is set to True when sending request message):

- **Middle range: from** (float64) is the upper limit of the lower range of the selected row
- **Lower range: value** (float64) is the value in the lower range of the selected row
- **Middle range: value** (float64) is the value in the middle range of the selected row
- **Middle range: to** (float64) is the lower limit of the upper range of the selected row
- **Upper range: value** (float64) is the value of the upper range of the selected row
- **Alternative range: value** (float64) is the value of the alternative range of the selected row
- **Error** described in the Response message>Body section

3DSwp.TimingEnable

Enables/disables the Advanced Timing in the 3D Sweeper module.

Arguments:

- **Enable** (unsigned int32) where 0=Disable, 1=Enable

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.TimingSend

Sends the Advanced Timing configuration of the 3D Sweeper module to the real time controller.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

3DSwp.FilePathsGet

Returns the list of file paths for the data saved by one single measurement (i.e. 1D and 2D sweeps save one single file, whereas a 3D sweep saves as many files as points configured for Step Channel 2).

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **File paths size** (int) is the size in bytes of the File paths string array
- **Number of file paths** (int) defines the number of elements of the File paths string array
- **File paths** (1D array string) returns the list of file paths. The size of each string item comes right before it as integer 32.
- **Error** described in the Response message>Body section

1D Sweeper

1DSwp.AcqChsSet

Arguments:

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of recorded channels. The indexes correspond to the list of Measurement in the Nanonis software.
To get the Measurements names use the *Signals.MeasNamesGet* function.
If this argument is empty, the software will take the Channel names as the valid argument to set the list of recorded channels.
- **Channel names size** (int) is the size in bytes of the channel names array. These are the names of the recorded channels
- **Channel names number** (int) is the number of elements of the channel names array
- **Channel names** (1D array string) is an array of channel names strings, where each string comes prepended by its size in bytes

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

1DSwp.AcqChsGet

Returns the list of recorded channels of the 1D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of the recorded channels. The indexes correspond to the list of Measurement in the Nanonis software.
- **Channels names size** (int) is the size in bytes of the Channels names string array
- **Number of channels** (int) is the number of elements of the Channels names string array
- **Channels names** (1D array string) returns the list of recorded channels names. The size of each string item comes right before it as integer 32
- **Error** described in the Response message>Body section

1DSwp.SwpSignalSet

Sets the Sweep signal in the 1D Sweeper.

Arguments:

- **Sweep channel name size** (int) is the number of characters of the sweep channel name string
- **Sweep channel name** (string) is the name of the signal selected for the sweep channel

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

1DSwp.SwpSignalGet

Returns the selected Sweep signal in the 1D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Sweep channel name size** (int) is the number of characters of the sweep channel name string
- **Sweep channel name** (string) is the name of the signal selected for the sweep channel
- **Error** described in the Response message>Body section

1DSwp.LimitsSet

Sets the limits of the Sweep signal in the 1D Sweeper.

Arguments:

- **Lower limit** (float32) defines the lower limit of the sweep range
- **Upper limit** (float32) defines the upper limit of the sweep range

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

1DSwp.LimitsGet

Returns the limits of the Sweep signal in the 1D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Lower limit** (float32) defines the lower limit of the sweep range
- **Upper limit** (float32) defines the upper limit of the sweep range
- **Error** described in the Response message>Body section

1DSwp.PropsSet

Sets the configuration of the parameters in the 1D Sweeper.

Arguments:

- **Initial Settling time (ms)** (float32)
- **Maximum slew rate (units/s)** (float32)
- **Number of steps** (int) defines the number of steps of the sweep. 0 points means no change
- **Period (ms)** (unsigned int16) where 0 means no change
- **Autosave** (int) defines if the sweep is automatically saved, where 0=no change, 1=On, 2=Off
- **Save dialog box** (int) defines if the save dialog box shows up or not, where 0=no change, 1=On, 2=Off
- **Settling time (ms)** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

1DSwp.PropsGet

Returns the configuration of the parameters in the 1D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Initial Settling time (ms)** (float32)
- **Maximum slew rate (units/s)** (float32)
- **Number of steps** (int) defines the number of steps of the sweep
- **Period (ms)** (unsigned int16)
- **Autosave** (unsigned int32) defines if the sweep is automatically saved, where 0=Off, 1=On
- **Save dialog box** (unsigned int32) defines if the save dialog box shows up or not, where 0=Off, 1=On
- **Settling time (ms)** (float32)
- **Error** described in the Response message>Body section

1DSwp.Start

Starts the sweep in the 1D Sweeper.

Arguments:

- **Get data** (unsigned int32) defines if the function returns the sweep data (1=True) or not (0=False)
- **Sweep direction** (unsigned int32) defines if the sweep starts from the lower limit (=1) or from the upper limit (=0)
- **Save base name string size** (int) defines the number of characters of the Save base name string
- **Save base name** (string) is the basename used by the saved files. If empty string, there is no change
- **Reset signal** (unsigned int32) where 0=Off, 1=On
- **Dummy value** (unsigned int16) used internally by the Nanonis software. It must be set to 0

Return arguments (if Send response back flag is set to True when sending request message):

- **Channels names size** (int) is the size in bytes of the Channels names string array
- **Number of channels** (int) is the number of elements of the Channels names string array
- **Channels names** (1D array string) returns the list of channels names. The size of each string item comes right before it as integer 32
- **Data rows** (int) defines the number of rows of the Data array
- **Data columns** (int) defines the number of columns of the Data array
- **Data** (2D array float32) returns the sweep data
- **Error** described in the Response message>Body section

1DSwp.Stop

Stops the sweep in the 1D Sweeper module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

1DSwp.Open

Opens the 1D Sweeper module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

1DSwp.ModeSet

Sets the sweep mode in the 1D Sweeper, which defines how the sweeper should react to each new data point set during the sweep.

Timed: Waits the amount of time set by the user (Settling Time) after setting a new sweep value using the Slew rate specified on the UI of the Sweeper. Does not explicitly wait for the target device to reach the value set.

Value Reached: Waits for feedback from the target device that the value set has actually been reached before starting integration. Does not require a Sweeper Slew Rate setting.

Continuous: After moving to start value, final value is set in one step and the measurement is performed continuously until the target has reached the specified value.

Arguments:

- **Sweep mode** (unsigned int16) where 0=Timed, 1=Value Reached, and 2=Continuous

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

1DSwp.ModeGet

Returns the sweep mode in the 1D Sweeper, which defines how the sweeper should react to each new data point set during the sweep.

Timed: Waits the amount of time set by the user (Settling Time) after setting a new sweep value using the Slew rate specified on the UI of the Sweeper. Does not explicitly wait for the target device to reach the value set.

Value Reached: Waits for feedback from the target device that the value set has actually been reached before starting integration. Does not require a Sweeper Slew Rate setting.

Continuous: After moving to start value, final value is set in one step and the measurement is performed continuously until the target has reached the specified value.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Sweep mode** (unsigned int16) where 0=Timed, 1=Value Reached, and 2=Continuous
- **Error** described in the Response message>Body section

1DSwp.SwpSignalListGet

Returns the list of names of available signals to sweep in the Generic Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Channels names size** (int) is the size in bytes of the Channels names string array
- **Number of channels** (int) is the number of elements of the Channels names string array
- **Channels names** (1D array string) returns the list of channels names. The size of each string item comes right before it as integer 32
- **Error** described in the Response message>Body section

High-Speed Sweeper

HSSwp.AcqChsSet

Sets the list of recorded channels of the High-Speed Sweeper.

Arguments:

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of recorded channels.
To obtain a list of the available channels, use the *HSSwp.AcqChsGet* function.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.AcqChsGet

Returns the list of recorded channels of the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of the recorded channels. The indexes correspond to the indices in the Available Channels indexes array.
- **Available Channels names size** (int) is the size in bytes of the available channels names array
- **Available Channels names number** (int) is the number of elements of the available channels names array
- **Available Channels names** (1D array string) returns an array of channel names strings, where each string comes prepended by its size in bytes
- **Number of available channels** (int) is the number of available channels. It defines the size of the Available Channels indexes array
- **Available Channels indexes** (1D array int) are the indexes of channels available for acquisition.
- **Error** described in the Response message>Body section

HSSwp.AutoReverseSet

Sets the auto-reverse configuration of the sweep axis in the High-Speed Sweeper.

Arguments:

- **On/Off** (int) defines if the auto-reverse functionality is on or off, where 0=Off, 1=On
- **Condition** (int) defines if the signal must be greater or less than the threshold for the reverse condition to activate, where 0 = >, 1 = <

- **Signal** (int) sets the signal for the reverse condition. The list of available signals is the same as the acquisition signals (see *HSSwp.AcqChsGet* function).
- **Threshold** (float32) defines the threshold to which the signal is compared.
- **Linkage to 1** (int) defines the linkage of the 2nd reverse condition to the first condition. Possible values: 0 = Off (no 2nd condition), 1 = OR (condition 1 or 2 must be met), 2 = AND (conditions 1 and 2 must be met at the same time), 3 = THEN (condition 1 must be met first, then condition 2).
- **Condition 2** (int) defines if the signal must be greater or less than the threshold for the 2nd reverse condition to activate, where 0 = >, 1 = <
- **Signal 2** (int) sets the signal for the 2nd reverse condition. The list of available signals is the same as the acquisition signals (see *HSSwp.AcqChsGet* function).
- **Threshold 2** (float32) defines the threshold to which the signal 2 is compared.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.AutoReverseGet

Returns the auto-reverse configuration of the sweep axis in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **On/Off** (int) specifies if the auto-reverse functionality is on or off, where 0=Off, 1=On
- **Condition** (int) specifies if the signal must be greater or less than the threshold for the reverse condition to activate, where 0 = >, 1 = <
- **Signal** (int) is the signal for the reverse condition. The list of available signals is the same as the acquisition signals (see *HSSwp.AcqChsGet* function).
- **Threshold** (float32) is the threshold to which the signal is compared.
- **Linkage to 1** (int) specifies the linkage of the 2nd reverse condition to the first condition. Possible values: 0 = Off (no 2nd condition), 1 = OR (condition 1 or 2 must be met), 2 = AND (conditions 1 and 2 must be met at the same time), 3 = THEN (condition 1 must be met first, then condition 2).
- **Condition 2** (int) specifies if the signal must be greater or less than the threshold for the 2nd reverse condition to activate, where 0 = >, 1 = <
- **Signal 2** (int) is the signal for the 2nd reverse condition. The list of available signals is the same as the acquisition signals (see *HSSwp.AcqChsGet* function).
- **Threshold 2** (float32) is the threshold to which the signal 2 is compared.
- **Error** described in the Response message>Body section

HSSwp.EndSettleSet

Sets the end settling time in the High-Speed Sweeper.

Arguments:

- **Threshold** (float32) defines the end settling time in seconds

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.EndSettlGet

Returns the end settling time in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Threshold** (float32) is the end settling time in seconds
- **Error** described in the Response message>Body section

HSSwp.NumSweepsSet

Sets the number of sweeps in the High-Speed Sweeper.

Arguments:

- **Number of sweeps** (unsigned int32) sets the number of sweeps (ignored when continuous is set)
- **Continuous** (int) sets the continuous sweep mode, where 0=Off, 1=On

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.NumSweepsGet

Returns the number of sweeps in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of sweeps** (unsigned int32) is the number of sweeps
- **Continuous** (int) specifies the continuous sweep mode, where 0=Off, 1=On
- **Error** described in the Response message>Body section

HSSwp.ResetSignalsSet

Specifies if the sweep and step signals should be reset to their initial values at the end of the sweep in the High-Speed Sweeper.

Arguments:

- **Reset Signals** (int) defines if the sweep and step signals are reset at the end of the sweep, where 0=Off, 1=On

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.ResetSignalsGet

Returns if the sweep and step signals are reset to their initial values at the end of the sweep in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Reset Signals** (int) returns if the sweep and step signals are reset at the end of the sweep, where 0=Off, 1=On
- **Error** described in the Response message>Body section

HSSwp.SaveBasenameSet

Sets the save basename in the High-Speed Sweeper.

Arguments:

- **Basename size** (int) is the size (number of characters) of the basename string
- **Basename** (string) is the base name used for the saved sweeps
- **Path size** (int) is the size (number of characters) of the path string
- **Path** (string) is the path used for the saved sweeps

Return arguments (if Send response back flag is set to True when sending request message to the server):

- **Error** described in the Response message>Body section

HSSwp.SaveBasenameGet

Returns the save basename in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message to the server):

- **Basename size** (int) is the size (number of characters) of the basename string
- **Basename** (string) is the base name used for the saved sweeps
- **Path size** (int) is the size (number of characters) of the path string
- **Path** (string) is the path used for the saved sweeps
- **Error** described in the Response message>Body section

HSSwp.SaveDataSet

Specifies if the data acquired in the High-Speed Sweeper is saved or not.

Arguments:

- **Save Data** (int) defines if the data is saved, where 0=Off, 1=On

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.SaveDataGet

Returns if the data acquired in the High-Speed Sweeper is saved or not.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Save Data** (int) returns if the data is saved, where 0=Off, 1=On
- **Error** described in the Response message>Body section

HSSwp.SaveOptionsSet

Sets save options in the High-Speed Sweeper.

Arguments:

- **Comment size** (int) is the size (number of characters) of the comment string
- **Comment** (string) is the comment saved in the header of the files. If empty string, there is no change
- **Modules names size** (int) is the size in bytes of the modules array. These are the modules whose parameters are saved in the header of the files
- **Modules names number** (int) is the number of elements of the modules names array
- **Modules names** (1D array string) is an array of modules names strings, where each string comes prepended by its size in bytes

HSSwp.SaveOptionsGet

Returns the saving options of the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message to the server):

- **Comment size** (int) is the size (number of characters) of the comment string
- **Comment** (string) is the comment saved in the header of the files
- **Modules parameters size** (int) is the size in bytes of the modules parameters array. These are the modules parameters saved in the header of the files
- **Modules parameters number** (int) is the number of elements of the modules parameters array
- **Modules parameters** (1D array string) is an array of modules names strings, where each string comes prepended by its size in bytes.
- **Error** described in the Response message>Body section

HSSwp.Start

Starts a sweep in the High-Speed Sweeper module.

When Send response back is set to True, it returns immediately afterwards.

Arguments:

- **Wait until done** (int) specifies whether the function waits with sending back the return arguments until the sweep is finished, where 0=Off (don't wait), 1=On (wait)
- **Timeout** (int) sets the wait timeout in milliseconds. Use -1 for indefinite wait. The Timeout is ignored when wait until done is off.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.Stop

Stops the sweep in the High-Speed Sweeper module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.StatusGet

Returns the status of the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Status** (unsigned int32) is status of the High-Speed Sweeper, where 0=Stopped, 1=Running
- **Error** described in the Response message>Body section

HSSwp.SwpChSigListGet

Returns the list of available signals for the Sweep channel of the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message to the server):

- **Signal names size** (int) is the size in bytes of the signal names array
- **Signal names number** (int) is the number of elements of the signal names array
- **Signal names** (1D array string) is an array of signal names strings, where each string comes prepended by its size in bytes.
- **Number of signals** (int) is the number of available sweep signals. It defines the size of the Available Channels indexes array
- **Signal indexes** (1D array int) are the indexes of signals available for the sweep channel.
- **Error** described in the Response message>Body section

HSSwp.SwpChSignalSet

Sets the Sweep Channel signal in the High-Speed Sweeper.

Arguments:

- **Sweep signal index** (int) is the index of the Sweep Signal. Use the *HSSwp.SwpChSigListGet* function to obtain a list of available sweep signals.
- **Timed Sweep** (int) enables or disables timed sweep mode. When on, the Sweep channel is ignored. 0=Off (sweep signal), 1=On (timed sweep)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.SwpChSignalGet

Returns the Sweep Channel signal in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Sweep signal index** (int) is the index of the Sweep Signal. Use the *HSSwp.SwpChSigListGet* function to obtain a list of available sweep signals.
- **Timed Sweep** (int) specifies if timed sweep mode is enabled, where 0=Off (sweep signal), 1=On (timed sweep)
- **Error** described in the Response message>Body section

HSSwp.SwpChLimitsSet

Sets the limits of the Sweep Channel in the High-Speed Sweeper.

Arguments:

- **Relative Limits** (int) specifies if the limits are absolute or relative to the current sweep signal value. Possible values are 0=Absolute limits, 1=Relative limits.
- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.SwpChLimitsGet

Returns the limits of the Sweep Channel in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Relative Limits** (int) specifies if the limits are absolute or relative to the current sweep signal value. Possible values are 0=Absolute limits, 1=Relative limits.
- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops
- **Error** described in the Response message>Body section

HSSwp.SwpChNumPtsSet

Sets the number of points for the Sweep Channel in the High-Speed Sweeper.

Arguments:

- **Number of points** (unsigned int32) sets the number of points of the sweep.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.SwpChNumPtsGet

Returns the number of points for the Sweep Channel in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of points** (int) returns the number of points of the sweep
- **Error** described in the Response message>Body section

HSSwp.SwpChTimingSet

Sets the timing parameters of the Sweep Channel in the High-Speed Sweeper.

Arguments:

- **Initial settling time (s)** (float32)
- **Settling time (s)** (float32)
- **Integration time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.SwpChTimingGet

Returns the timing parameters of the Sweep Channel in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Initial settling time (s)** (float32)
- **Settling time (s)** (float32)
- **Integration time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)
- **Error** described in the Response message>Body section

HSSwp.SwpChBwdSwSet

Enables or disables the backward sweep for the sweep channel in the High-Speed Sweeper.

Arguments:

- **Bwd Sweep** (unsigned int32) defines if the backward sweep is enabled, where 0=Off, 1=On

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.SwpChBwdSwGet

Returns if the backward sweep of the sweep channel in the High-Speed Sweeper is enabled or not.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Bwd Sweep** (unsigned int32) specifies if the backward sweep is enabled, where 0=Off, 1=On
- **Error** described in the Response message>Body section

HSSwp.SwpChBwdDelaySet

Sets the delay between forward and backward sweep of the sweep channel in the High-Speed Sweeper.

Arguments:

- **Bwd Delay** (float32) sets the delay between forward and backward sweep in seconds.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.SwpChBwdDelayGet

Returns the delay between forward and backward sweep of the sweep channel in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Bwd Delay** (float32) is the delay between forward and backward sweep in seconds
- **Error** described in the Response message>Body section

HSSwp.StpChSigListGet

Returns the list of available signals for the Step channel of the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message to the server):

- **Signal names size** (int) is the size in bytes of the signal names array
- **Signal names number** (int) is the number of elements of the signal names array
- **Signal names** (1D array string) is an array of signal names strings, where each string comes prepended by its size in bytes.
- **Number of signals** (int) is the number of available sweep signals. It defines the size of the Available Channels indexes array
- **Signal indexes** (1D array int) are the indexes of signals available for the sweep channel.
- **Error** described in the Response message>Body section

HSSwp.StpChSignalSet

Sets the Step Channel signal in the High-Speed Sweeper.

Arguments:

- **Step signal index** (int) is the index of the Step Signal. Use the *HSSwp.StpChSigListGet* function to obtain a list of available step signals.
- **Unused** (int) enables or disables the step channel, where 0=Off (step channel disabled), 1=On (step channel enabled)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.StpChSignalGet

Returns the Step Channel signal in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Step signal index** (int) is the index of the Step Signal. Use the *HSSwp.StpChSigListGet* function to obtain a list of available step signals.
- **Unused** (int) specifies if the step channel is enabled, where 0=Off (step channel disabled), 1=On (step channel enabled)
- **Error** described in the Response message>Body section

HSSwp.StpChLimitsSet

Sets the limits of the Step Channel in the High-Speed Sweeper.

Arguments:

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.StpChLimitsGet

Returns the limits of the Step Channel in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Start** (float32) defines the value where the sweep starts
- **Stop** (float32) defines the value where the sweep stops
- **Error** described in the Response message>Body section

HSSwp.StpChNumPtsSet

Sets the number of points for the Step Channel in the High-Speed Sweeper.

Arguments:

- **Number of points** (unsigned int32) sets the number of points of the sweep.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.StpChNumPtsGet

Returns the number of points for the Step Channel in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of points** (int) returns the number of points of the sweep
- **Error** described in the Response message>Body section

HSSwp.StpChTimingSet

Sets the timing parameters of the Step Channel in the High-Speed Sweeper.

Arguments:

- **Initial settling time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

HSSwp.StpChTimingGet

Returns the timing parameters of the Step Channel in the High-Speed Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Initial settling time (s)** (float32)
- **Maximum slew rate (units/s)** (float32)
- **Error** described in the Response message>Body section

Signals

Signals.NamesGet

Returns the signals names list of the 128 signals available in the software.

The 128 signals are physical inputs, physical outputs and internal channels. By searching in the list the channel's name you are interested in, you can get its index (0-127).

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Signals names size** (int) is the size in bytes of the signals names array
- **Signals names number** (int) is the number of elements of the signals names array
- **Signals names** (1D array string) returns an array of signals names strings, where each string comes prepended by its size in bytes
- **Error** described in the Response message>Body section

Signals.CalibrGet

Returns the calibration and offset of the selected signal.

Arguments:

- **Signal index** (int) is comprised between 0 and 127

Return arguments (if Send response back flag is set to True when sending request message):

- **Calibration per volt** (float32)
- **Offset in physical units** (float32)
- **Error** described in the Response message>Body section

Signals.RangeGet

Returns the range limits of the selected signal.

Arguments:

- **Signal index** (int) is comprised between 0 and 127

Return arguments (if Send response back flag is set to True when sending request message):

- **Maximum limit** (float32)
- **Minimum limit** (float32)
- **Error** described in the Response message>Body section

Signals.ValGet

Returns the current value of the selected signal (oversampled during the Acquisition Period time, Tap).

Signal measurement principle:

The signal is continuously oversampled with the Acquisition Period time, Tap, specified in the TCP receiver module. Every Tap second, the oversampled data is "published". This VI function waits for the next oversampled data to be published and returns its value. Calling this function does not trigger a signal measurement; it waits for data to be published! Thus, this function returns a value 0 to Tap second after being called.

An important consequence is that if you change a signal and immediately call this function to read a measurement you might get "old" data (i.e. signal data measured before you changed the signal). The solution to get only new data is to set Wait for newest data to True. In this case, the first published data is discarded and only the second one is returned.

Arguments:

- **Signal index** (int) is comprised between 0 and 127
- **Wait for newest data** (unsigned int32) selects whether the function returns the next available signal value or if it waits for a full period of new data. If False, this function returns a value 0 to Tap seconds after being called. If True, the function discard the first oversampled signal value received but returns the second value received. Thus, the function returns a value Tap to 2*Tap seconds after being called. It could be 0=False or 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Signal value** (float32) is the value of the selected signal in physical units
- **Error** described in the Response message>Body section

Signals.ValsGet

Returns the current values of the selected signals (oversampled during the Acquisition Period time, Tap).

Arguments:

- **Signals indexes size** (int) is the size of the Signals indexes array
- **Signals indexes** (1D array int) sets the selection of signals indexes, comprised between 0 and 127
- **Wait for newest data** (unsigned int32) selects whether the function returns the next available signal value or if it waits for a full period of new data. If False, this function returns a value 0 to Tap seconds after being called. If True, the function discard the first oversampled signal value received but returns the second value received. Thus, the function returns a value Tap to 2*Tap seconds after being called. It could be 0=False or 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Signals values size** (int) is the size of the Signals values array
- **Signals values** (1D array float32) returns the values of the selected signals in physical units
- **Error** described in the Response message>Body section

Signals.MeasNamesGet

Returns the list of measurement channels names available in the software.

Important Note: The Measurement channels don't correspond to the Signals. Measurement channels are used in sweepers whereas the Signals are used by the graphs and other modules.

By searching in the list the channels's names you are interested in, you can know its index. This index is then used e.g. to get/set the recorded channels in Sweepers, for example by using the *GenSwp.ChannelsGet* and *GenSwp.ChannelsSet* functions for the 1D Sweeper.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Measurement channels list size** (int) is the size in bytes of the Measurement channels list array
- **Number of Measurement channels** (int) is the number of elements of the Measurement channels list array
- **Measurement channels list** (1D array string) returns an array of names, where each array element is preceded by its size in bytes
- **Error** described in the Response message>Body section

Signals.AddRTGet

Returns the list of names of additional RT signals available, and the names of the signals currently assigned to the Internal 23 and 24 signals.

This can be found in the Signals Manager.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Additional RT signals names size** (int) is the size in bytes of the Additional RT signals names array
- **Number of Additional RT signals** (int) is the number of elements of the Additional RT signals names array
- **Additional RT signals names** (1D array string) returns the list of additional RT signals which can be assigned to Internal 23 and 24. Each array element is preceded by its size in bytes
- **Additional RT signal 1** (string) is the name of the RT signal assigned to the Internal 23 signal
- **Additional RT signal 2** (string) is the name of the RT signal assigned to the Internal 24 signal
- **Error** described in the Response message>Body section

Signals.AddRTSet

Assigns additional RT signals to the Internal 23 and 24 signals in the Signals Manager.

This function links advanced RT signals to Internal 23 and Internal 24.

Arguments:

- **Additional RT signal 1** (int) is the index of the RT signal assigned to the Internal 23 signal
- **Additional RT signal 2** (int) is the index of the RT signal assigned to the Internal 24 signal

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

User Inputs

UserIn.CalibrSet

Sets the calibration of the selected user input.

Arguments:

- **Input index** (int) sets the input to be used, where index could be any value from 1 to the available inputs
- **Calibration per volt** (float32)
- **Offset in physical units** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

User Outputs

UserOut.ModeSet

Sets the mode (User Output, Monitor, Calculated signal) of the selected user output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs
- **Output mode** (unsigned int16) sets the output mode of the selected output, where 0=User Output, 1=Monitor, 2=Calc.Signal

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

UserOut.ModeGet

Returns the mode (User Output, Monitor, Calculated signal) of the selected user output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs

Return arguments (if Send response back flag is set to True when sending request message):

- **Output mode** (unsigned int16) returns the output mode of the selected output, where 0=User Output, 1=Monitor, 2=Calc.Signal, 3=Override
- **Error** described in the Response message>Body section

UserOut.MonitorChSet

Sets the monitor channel of the selected output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs
- **Monitor channel index** (int) sets the index of the channel to monitor. The index is comprised between 0 and 127 for the physical inputs, physical outputs, and internal channels. To see which signal has which index, see *Signals.NamesGet* function

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

UserOut.MonitorChGet

Returns the monitor channel of the selected output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs

Return arguments (if Send response back flag is set to True when sending request message):

- **Monitor channel index** (int) returns the index of the channel to monitor. The index is comprised between 0 and 127 for the physical inputs, physical outputs, and internal channels. To see which signal has which index, see *Signals.NamesGet* function
- **Error** described in the Response message>Body section

UserOut.ValSet

Sets the value of the selected user output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs
- **Output value** (float32) is the value applied to the selected user output in physical units

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

UserOut.CalibrSet

Sets the calibration of the selected user output or monitor channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs
- **Calibration per volt** (float32)
- **Offset in physical units** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

UserOut.CalcSignalNameSet

Sets the Calculated Signal name of the selected output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs
- **Calculated signal name size** (int) is the number of characters of the Calculated signal name string
- **Calculated signal name** (string) is the name of the calculated signal configured for the selected output

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

UserOut.CalcSignalNameGet

Returns the Calculated Signal name of the selected output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs

Return arguments (if Send response back flag is set to True when sending request message):

- **Calculated signal name size** (int) is the number of characters of the Calculated signal name string
- **Calculated signal name** (string) is the name of the calculated signal configured for the selected output
- **Error** described in the Response message>Body section

UserOut.CalcSignalConfigSet

Sets the configuration of the Calculated Signal for the selected output channel.

The configuration is a combination of 4 parts that creates a formula. Each part of the formula is a parameter/math operation and a value (which depending on the parameter/math operation is applicable or not).

The possible values for the parameter/math operation of part 1 are:

0=None, 5=Constant, 10=Signal Index

The possible values for the parameter/math operation of parts 2, 3 and 4 are:

0=None, 1=Add Constant, 1=Subtract Constant, 3=Multiply Constant, 4=Divide Constant, 6=Add Signal, 7=Subtract Signal, 8=Multiply Signal, 9=Divide Signal, 11=Exponent, 12=Absolute, 13= Negate, 14= Log

There is no mathematical operator precedence in operation here; the equations are executed in a strict left-to-right (part 1 to part 4) fashion.

This is equivalent to defining the calculations as (((Part 1) Part 2) Part 3) Part 4).

For example:

The average of Input 1 and Input 2 is defined as “Input 1 + Input 2 / 2”.

The sum of Input 1 plus half of Input 2 is defined as “Input 2 / 2 + Input 1”.

The reciprocal of Input 1 is defined as “1 / Input 1”.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs
- **Operation 1** (unsigned int16) is the parameter or math operation selected as the 1st part of the configuration formula.
- **Value 1** (float32) is a constant value or signal index, depending on the precedent operation
- **Operation 2** (unsigned int16) is the parameter or math operation selected as the 2nd part of the configuration formula.
- **Value 2** (float32) is a constant value or signal index, depending on the precedent operation
- **Operation 3** (unsigned int16) is the parameter or math operation selected as the 3rd part of the configuration formula.
- **Value 3** (float32) is a constant value or signal index, depending on the precedent operation
- **Operation 4** (unsigned int16) is the parameter or math operation selected as the 4th part of the configuration formula.
- **Value 4** (float32) is a constant value or signal index, depending on the precedent operation

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

UserOut.CalcSignalConfigGet

Returns the configuration of the Calculated Signal for the selected output channel.

The configuration is a combination of 4 parts that creates a formula. Each part of the formula is a parameter/math operation and a value (which depending on the parameter/math operation is applicable or not).

The possible values for the parameter/math operation of part 1 are:

0=None, 5=Constant, 10=Signal Index

The possible values for the parameter/math operation of parts 2, 3 and 4 are:

0=None, 1=Add Constant, 1=Subtract Constant, 3=Multiply Constant, 4=Divide Constant, 6=Add Signal, 7=Subtract Signal, 8=Multiply Signal, 9=Divide Signal, 11=Exponent, 12=Absolute, 13= Negate, 14= Log

There is no mathematical operator precedence in operation here; the equations are executed in a strict left-to-right (part 1 to part 4) fashion.

This is equivalent to defining the calculations as (((Part 1) Part 2) Part 3) Part 4).

For example:

The average of Input 1 and Input 2 is defined as “Input 1 + Input 2 / 2”.

The sum of Input 1 plus half of Input 2 is defined as “Input 2 / 2 + Input 1”.

The reciprocal of Input 1 is defined as “1 / Input 1”.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs

Return arguments (if Send response back flag is set to True when sending request message):

- **Operation 1** (unsigned int16) is the parameter or math operation selected as the 1st part of the configuration formula.
- **Value 1** (float32) is a constant value or signal index, depending on the precedent operation
- **Operation 2** (unsigned int16) is the parameter or math operation selected as the 2nd part of the configuration formula.
- **Value 2** (float32) is a constant value or signal index, depending on the precedent operation
- **Operation 3** (unsigned int16) is the parameter or math operation selected as the 3rd part of the configuration formula.
- **Value 3** (float32) is a constant value or signal index, depending on the precedent operation
- **Operation 4** (unsigned int16) is the parameter or math operation selected as the 4th part of the configuration formula.
- **Value 4** (float32) is a constant value or signal index, depending on the precedent operation
- **Error** described in the Response message>Body section

UserOut.LimitsSet

Sets the physical limits (in calibrated units) of the selected output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs
- **Upper limit** (float32) defines the upper physical limit of the user output
- **Lower limit** (float32) defines the lower physical limit of the user output
- **Raw limits?** (unsigned int32) selects whether to set the physical limits (0) or the raw limits (1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

UserOut.LimitsGet

Returns the physical limits (in calibrated units) of the selected output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs
- **Raw limits?** (unsigned int32) selects whether to get the physical limits (0) or the raw limits (1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Upper limit** (float32) defines the upper physical limit of the user output
- **Lower limit** (float32) defines the lower physical limit of the user output
- **Error** described in the Response message>Body section

UserOut.SlewRateSet

Sets the slew rate (in calibrated units) of the selected output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs
- **Slew Rate** (float64) defines the calibrated slew rate of the user output

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

UserOut.SlewRateGet

Returns the slew rate (in calibrated units) of the selected output channel.

Arguments:

- **Output index** (int) sets the output to be used, where index could be any value from 1 to the number of available outputs

Return arguments (if Send response back flag is set to True when sending request message):

- **Slew Rate** (float64) is the calibrated slew rate of the user output
- **Error** described in the Response message>Body section

Digital Lines

DigLines.PropsSet

Configures the properties of a digital line.

Arguments:

- **Digital line** (unsigned int32) defines the line to configure, from 1 to 8
- **Port** (unsigned int32) selects the digital port, where 0=Port A, 1=Port B, 2=Port C, 3=Port D, 4=Expanded Port 0, 5=Expanded Port 1, 6=Expanded Port 2, etc
The expanded DIO ports are available on selected systems through a dedicated DIO card.
- **Direction** (unsigned int32) is the direction of the selected digital line, where 0=Input, 1=Output
- **Polarity** (unsigned int32), where 0=Low active, 1=High active

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

DigLines.OutStatusSet

Sets the status of a digital output line.

Arguments:

- **Port** (unsigned int32) selects the digital port, where 0=Port A, 1=Port B, 2=Port C, 3=Port D, 4=Expanded Port 0, 5=Expanded Port 1, 6=Expanded Port 2, etc
The expanded DIO ports are available on selected systems through a dedicated DIO card.
- **Digital line** (unsigned int32) defines the output line to configure, from 1 to 8
- **Status** (unsigned int32) sets whether the output is active or inactive, where 0=Inactive, 1=Active

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

DigLines.TTLValGet

Reads the actual TTL voltages present at the pins of the selected port.

Arguments:

- **Port** (unsigned int32) selects the digital port, where 0=Port A, 1=Port B, 2=Port C, 3=Port D, 4=Expanded Port 0, 5=Expanded Port 1, 6=Expanded Port 2, etc
The expanded DIO ports are available on selected systems through a dedicated DIO card.

Return arguments (if Send response back flag is set to True when sending request message):

- **TTL voltages size** (int) is the size of the TTL voltages array
- **TTL voltages** (1D array unsigned int32) sets whether the output is active or inactive, where 0=Inactive, 1=Active
- **Error** described in the Response message>Body section

DigLines.Pulse

Configures and starts the pulse generator on the selected digital outputs.

Arguments:

- **Port** (unsigned int32) selects the digital port, where 0=Port A, 1=Port B, 2=Port C, 3=Port D, 4=Expanded Port 0, 5=Expanded Port 1, 6=Expanded Port 2, etc
The expanded DIO ports are available on selected systems through a dedicated DIO card.
- **Digital lines size** (int) is the size of the Digital lines array
- **Digital lines** (1D array unsigned int8) defines the output lines to pulse, from 1 to 8
- **Pulse width (s)** (float32) defines how long the outputs are active
- **Pulse pause (s)** (float32) defines how long the outputs are inactive
- **Number of pulses** (int) defines how many pulses to generate, where valid values are from 1 to 32767
- **Wait until finished** (unsigned int32), if True this function waits until all pulses have been generated before the response is sent back, where 0=False, 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Data Logger

DataLog.Open

Opens the Data Logger module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

DataLog.Start

Starts the acquisition in the Data Logger module.

Before using this function, select the channels to record in the Data Logger.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

DataLog.Stop

Stops the acquisition in the Data Logger module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

DataLog.StatusGet

Returns the status parameters from the Data Logger module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Start time size** (int) returns the number of bytes corresponding to the Start time string
- **Start time** (string) returns a timestamp of the moment when the acquisition started
- **Acquisition elapsed hours** (unsigned int16) returns the number of hours already passed since the acquisition started
- **Acquisition elapsed minutes** (unsigned int16) returns the number of minutes displayed on the Data Logger
- **Acquisition elapsed seconds** (float32) returns the number of seconds displayed on the Data Logger
- **Stop time size** (int) returns the number of bytes corresponding to the Stop time string
- **Stop time** (string) returns a timestamp of the moment when the acquisition Stopped
- **Saved file path size** (int) returns the number of bytes corresponding to the Saved file path string
- **Saved file path** (string) returns the path of the last saved file
- **Points counter** (int) returns the number of points (averaged samples) to save into file.
This parameter updates while running the acquisition
- **Error** described in the Response message>Body section

DataLog.ChsSet

Sets the list of recorded channels in the Data Logger module.

Arguments:

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of recorded channels. The index are comprised between 0 and 127, and it corresponds to the full list of signals available in the system.
To get the signal name and its corresponding index in the list of the 128 available signals in the Nanonis Controller, use the *Signal.NamesGet* function, or check the RT Idx value in the Signals Manager module

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

DataLog.ChsGet

Returns the list of recorded channels in the Data Logger module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of recorded channels. The indexes are comprised between 0 and 127, and it corresponds to the full list of signals available in the system.
To get the signal name and its corresponding index in the list of the 128 available signals in the Nanonis Controller, use the *Signal.NamesGet* function, or check the RT Idx value in the Signals Manager module
- **Error** described in the Response message>Body section

DataLog.PropsSet

Sets the acquisition configuration and the save options in the Data Logger module.

Arguments:

- **Acquisition mode** (unsigned int16) means that if Timed (=2), the selected channels are acquired during the acquisition duration time or until the user presses the Stop button.
If Continuous (=1), the selected channels are acquired continuously until the user presses the Stop button.
If 0, there is no change in the acquisition mode.
The acquired data are saved every time the averaged samples buffer reaches 25.000 samples and when the acquisition stops
- **Acquisition duration(hours)** (int) sets the number of hours the acquisition should last. Value -1 means no change
- **Acquisition duration (minutes)** (int) sets the number of minutes. Value -1 means no change
- **Acquisition duration (seconds)** (float32) sets the number of seconds. Value -1 means no change
- **Averaging** (int) sets how many data samples (received from the real-time system) are averaged for one data point saved into file. By increasing this value, the noise might decrease, and fewer points per seconds are recorded.
Use 0 to skip changing this parameter
- **Basename size** (int) is the size in bytes of the Basename string
- **Basename** (string) is base name used for the saved images
- **Comment size** (int) is the size in bytes of the Comment string
- **Comment** (string) is comment saved in the file
- **Size of the list of modules** (int) is the size in bytes of the List of modules string array
- **Number of modules** (int) is the number of elements of the List of modules string array
- **List of modules** (1D array string) sets the modules names whose parameters will be saved in the header of the files. The size of each string item should come right before it as integer 32

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

DataLog.PropsGet

Returns the acquisition configuration and the save options in the Data Logger module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Acquisition mode** (unsigned int16) means that if Timed (=1), the selected channels are acquired during the acquisition duration time or until the user presses the Stop button.
If Continuous (=0), the selected channels are acquired continuously until the user presses the Stop button.
The acquired data are saved every time the averaged samples buffer reaches 25.000 samples and when the acquisition stops
- **Acquisition duration(hours)** (int) returns the number of hours the acquisition lasts
- **Acquisition duration (minutes)** (int) returns the number of minutes
- **Acquisition duration (seconds)** (float32) returns the number of seconds
- **Averaging** (int) returns how many data samples (received from the real-time system) are averaged for one data point saved into file
- **Basename size** (int) returns the size in bytes of the Basename string
- **Basename** (string) returns the base name used for the saved images
- **Comment size** (int) returns the size in bytes of the Comment string
- **Comment** (string) returns the comment saved in the file
- **Error** described in the Response message>Body section

TCP Logger

TCPLog.Start

Starts the acquisition in the TCP Logger module.

Before using this function, select the channels to record in the TCP Logger.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

TCPLog.Stop

Stops the acquisition in the TCP Logger module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

TCPLog.ChsSet

Sets the list of recorded channels in the TCP Logger module.

Arguments:

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of recorded channels. The indexes are comprised between 0 and 127, and it corresponds to the full list of signals available in the system.
To get the signal name and its corresponding index in the list of the 128 available signals in the Nanonis Controller, use the *Signal.NamesGet* function, or check the RT Idx value in the Signals Manager module

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

TCPLog.OversampleSet

Sets the oversampling value in the TCP Logger.

Arguments:

- **Oversampling value** (int) sets the oversampling index, where index could be any value from 0 to 1000

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

TCPLog.StatusGet

Returns the current status of the TCP Logger.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Status** (int) returns an index which corresponds to one of the following status: 0=disconnected, 1=idle, 2=start, 3=stop, 4=running, 5=TCP connect, 6=TCP disconnect, 7=buffer overflow
- **Error** described in the Response message>Body section

Oscilloscope High Resolution

OsciHR.ChSet

Sets the measured signal index of the selected channel from the Oscilloscope High Resolution.

Arguments:

- **Channel index** (int) sets the oscilloscope channel to be used, where index could be any value from 1 to 4 when the Oscilloscope supports 4 channels. Otherwise the channel index should be set to 0 for the 1-channel version of this graph
- **Signal index** (int) sets the signal to be used, where index could be any value from 0 to 15

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.ChGet

Returns the measured signal index of the selected channel from the Oscilloscope High Resolution.

Arguments:

- **Channel index** (int) sets the oscilloscope channel to be used, where index could be any value from 1 to 4 when the Oscilloscope supports 4 channels. Otherwise the channel index should be set to 0 for the 1-channel version of this graph

Return arguments (if Send response back flag is set to True when sending request message):

- **Signal index** (int) returns the signal index
- **Error** described in the Response message>Body section

OsciHR.OversampleSet

Sets the oversampling index of the Oscilloscope High Resolution.

Choosing to acquire data at lower rate than the maximum 1MS/s allows for an improved S/N ratio and also increases the time window for the acquisition for a given number of samples.

Arguments:

- **Oversampling index** (int) sets the oversampling index, where index could be any value from 0 to 10

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.OversampleGet

Returns the oversampling index of the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Oversampling index** (int) gets the oversampling index, where index could be any value from 0 to 10
- **Error** described in the Response message>Body section

OsciHR.CalibrModeSet

Sets the calibration mode of the selected channel from the Oscilloscope High Resolution.

Select between Raw Values or Calibrated Values. This setting affects the data displayed in the graph, and trigger level and hysteresis values.

Arguments:

- **Channel index** (int) sets the oscilloscope channel to be used, where index could be any value from 1 to 4 when the Oscilloscope supports 4 channels. Otherwise the channel index should be set to 0 for the 1-channel version of this graph
- **Calibration mode** (unsigned int16), where 0=Raw values and 1=Calibrated values

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.CalibrModeGet

Returns the calibration mode of the selected channel from the Oscilloscope High Resolution.

Arguments:

- **Channel index** (int) sets the oscilloscope channel to be used, where index could be any value from 1 to 4 when the Oscilloscope supports 4 channels. Otherwise the channel index should be set to 0 for the 1-channel version of this graph

Return arguments (if Send response back flag is set to True when sending request message):

- **Calibration mode** (unsigned int16), where 0=Raw values and 1=Calibrated values
- **Error** described in the Response message>Body section

OsciHR.SamplesSet

Sets the number of samples to acquire in the Oscilloscope High Resolution.

Arguments:

- **Number of samples** (int)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.SamplesGet

Returns the number of samples to acquire in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of samples** (int)
- **Error** described in the Response message>Body section

OsciHR.PreTrigSet

Sets the Pre-Trigger Samples or Seconds in the Oscilloscope High Resolution.

If Pre-Trigger (s) is NaN or Inf or below 0, Pre-Trigger Samples is taken into account instead of seconds.

Arguments:

- **Pre-Trigger samples** (unsigned int32)
- **Pre-Trigger (s)** (float64)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.PreTrigGet

Returns the Pre-Trigger Samples in the Oscilloscope High Resolution.

If Pre-Trigger (s) is NaN or Inf or below 0, Pre-Trigger Samples is taken into account instead of seconds.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Pre-Trigger samples** (int)
- **Error** described in the Response message>Body section

OsciHR.Run

Starts the Oscilloscope High Resolution module.

The Oscilloscope High Resolution module does not run when its front panel is closed. To automate measurements it might be required to run the module first using this function.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.OsciDataGet

Returns the graph data of the selected channel from the Oscilloscope High Resolution.

Arguments:

- **Channel index** (int) sets the oscilloscope channel to be used, where index could be any value from 1 to 4 when the Oscilloscope supports 4 channels. Otherwise the channel index should be set to 0 for the 1-channel version of this graph
- **Data to get** (unsigned int16), where 0=Current returns the currently displayed data and 1=Next trigger waits for the next trigger to retrieve data
- **Timeout (s)** (float64), where -1 means waiting forever

Return arguments (if Send response back flag is set to True when sending request message):

- **Data t0 size** (int) is the number of characters of the Data t0 string
- **Data t0** (string) is the timestamp of the 1st acquired point
- **Data dt** (float64) is the time distance between two acquired points
- **Data Y size** (int) is the number of data points in Data Y
- **Data Y** (1D array float32) is the data acquired in the oscilloscope
- **Timeout** (unsigned int32) is 0 when no timeout occurred, and 1 when a timeout occurred
- **Error** described in the Response message>Body section

OsciHR.TrigModeSet

Sets the trigger mode in the Oscilloscope High Resolution.

Arguments:

- **Trigger mode** (unsigned int16), 0=Immediate means triggering immediately whenever the current data set is received by the host software, 1=Level where the trigger detection is performed on the non-averaged raw channel data (1MS/s), and 2=Digital where the trigger detection on the LS-DIO channels is performed at 500kS/s. Trigger detection on the HS-DIO channels is performed at 10MS/s

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.TrigModeGet

Returns the trigger mode in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Trigger mode** (unsigned int16), where 0=Immediate, 1=Level, and 2=Digital
- **Error** described in the Response message>Body section

OsciHR.TrigLevChSet

Sets the Level Trigger Channel index in the Oscilloscope High Resolution.

Trigger detection is performed on the non-averaged raw channel data.

Arguments:

- **Level trigger channel index** (int)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.TrigLevChGet

Returns the Level Trigger Channel index in the Oscilloscope High Resolution.

Trigger detection is performed on the non-averaged raw channel data.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Level trigger channel index** (int)
- **Error** described in the Response message>Body section

OsciHR.TrigLevValSet

Sets the Level Trigger value in the Oscilloscope High Resolution.

Arguments:

- **Level trigger value** (float64)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.TrigLevValGet

Returns the Level Trigger value in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Level trigger value** (float64)
- **Error** described in the Response message>Body section

OsciHR.TrigLevHystSet

Sets the Level Trigger Hysteresis in the Oscilloscope High Resolution.

Arguments:

- **Level trigger Hysteresis** (float64)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.TrigLevHystGet

Returns the Level Trigger Hysteresis in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Level trigger Hysteresis** (float64)
- **Error** described in the Response message>Body section

OsciHR.TrigLevSlopeSet

Sets the Level Trigger Slope in the Oscilloscope High Resolution.

Arguments:

- **Level trigger slope** (unsigned int16), where 0=Rising and 1=Falling

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.TrigLevSlopeGet

Returns the Level Trigger Slope in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Level trigger slope** (unsigned int16), where 0=Rising and 1=Falling
- **Error** described in the Response message>Body section

OsciHR.TrigDigChSet

Sets the Digital Trigger Channel in the Oscilloscope High Resolution.

Trigger detection on the LS-DIO channels is performed at 500kS/s. Trigger detection on the HS-DIO channels is performed at 10MS/s.

Arguments:

- **Digital trigger channel index** (int), where index can be any value from 0 to 35. Low Speed Port A lines are indexes 0 to 7, Low Speed Port B lines are indexes 8 to 15, Low Speed Port C lines are indexes 16 to 23, Low Speed Port D lines are indexes 24 to 31, and High Speed Port lines are indexes 32 to 35

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.TrigDigChGet

Returns the Digital Trigger Channel in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Digital trigger channel index** (int), where index can be any value from 0 to 35. Low Speed Port A lines are indexes 0 to 7, Low Speed Port B lines are indexes 8 to 15, Low Speed Port C lines are indexes 16 to 23, Low Speed Port D lines are indexes 24 to 31, and High Speed Port lines are indexes 32 to 35
- **Error** described in the Response message>Body section

OsciHR.TrigArmModeSet

Sets the Trigger Arming Mode in the Oscilloscope High Resolution.

Arguments:

- **Trigger arming mode** (unsigned int16), where 0=Single shot means recording the next available data and stopping acquisition. and 1=Continuous means recording every available data and automatically re-triggers the acquisition

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.TrigArmModeGet

Returns the Trigger Arming Mode in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Trigger arming mode** (unsigned int16), where 0=Single shot means recording the next available data and stopping acquisition. and 1=Continuous means recording every available data and automatically re-triggers the acquisition
- **Error** described in the Response message>Body section

OsciHR.TrigDigSlopeSet

Sets the Digital Trigger Slope in the Oscilloscope High Resolution.

Arguments:

- **Digital trigger slope** (unsigned int16), where 0=Rising and 1=Falling

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.TrigDigSlopeGet

Returns the Digital Trigger Slope in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Digital trigger slope** (unsigned int16), where 0=Rising and 1=Falling
- **Error** described in the Response message>Body section

OsciHR.TrigRearm

Rearms the trigger in the Oscilloscope High Resolution module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.PSDShow

Shows or hides the PSD section of the Oscilloscope High Resolution.

Arguments:

- **Show PSD section** (unsigned int32), where 0=Hide and 1=Show

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.PSDWeightSet

Sets the PSD Weighting in the Oscilloscope High Resolution.

Arguments:

- **PSD Weighting** (unsigned int16), where 0=Linear means that the averaging combines Count spectral records with equal weighting and then stops, whereas 1=Exponential means that the averaging process is continuous and new spectral data have a higher weighting than older ones

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.PSDWeightGet

Returns the PSD Weighting in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **PSD Weighting** (unsigned int16), where 0=Linear means that the averaging combines Count spectral records with equal weighting and then stops, whereas 1=Exponential means that the averaging process is continuous and new spectral data have a higher weighting than older ones
- **Error** described in the Response message>Body section

OsciHR.PSDWindowSet

Sets the PSD Window Type in the Oscilloscope High Resolution.

Arguments:

- **PSD window type** (unsigned int16) is the window function applied to the timed signal before calculating the power spectral density, where 0=None, 1=Hanning, 2=Hamming, etc

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.PSDWindowGet

Returns the PSD Window Type in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **PSD window type** (unsigned int16) is the window function applied to the timed signal before calculating the power spectral density, where 0=None, 1=Hanning, 2=Hamming, etc
- **Error** described in the Response message>Body section

OsciHR.PSDAvgTypeSet

Sets the PSD Averaging Type in the Oscilloscope High Resolution.

Arguments:

- **PSD averaging type** (unsigned int16), where 0=None, 1=Vector, 2=RMS, 3=Peak hold

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.PSDAvgTypeGet

Returns the PSD Averaging Type in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **PSD averaging type** (unsigned int16), where 0=None, 1=Vector, 2=RMS, 3=Peak hold
- **Error** described in the Response message>Body section

OsciHR.PSDAavgCountSet

Sets the PSD Averaging Count used by the RMS and Vector averaging types in the Oscilloscope High Resolution.

Arguments:

- **PSD averaging count** (int)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.PSDAavgCountGet

Returns the PSD Averaging Count used by the RMS and Vector averaging types in the Oscilloscope High Resolution.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **PSD averaging count** (int)
- **Error** described in the Response message>Body section

OsciHR.PSDAavgRestart

Restarts the PSD averaging process in the Oscilloscope High Resolution module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

OsciHR.PSDDataGet

Returns the Power Spectral Density data from the Oscilloscope High Resolution.

Arguments:

- **Data to get** (unsigned int16), where 0=Current returns the currently displayed data and 1=Next trigger waits for the next trigger to retrieve data
- **Timeout (s)** (float64), where -1 means waiting forever

Return arguments (if Send response back flag is set to True when sending request message):

- **Data f0** (float64) is the x coordinate of the 1st acquired point
- **Data df** (float64) is the frequency distance between two acquired points
- **Data Y size** (int) is the number of data points in Data Y
- **Data Y** (1D array float64) is the PSD data acquired in the oscilloscope
- **Timeout** (unsigned int32) is 0 when no timeout occurred, and 1 when a timeout occurred
- **Error** described in the Response message>Body section

Script

Script.Open

Opens the Script module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.Load

Loads a script in the script module.

Arguments:

- **Script index** (int) sets the slot where the script will be loaded and covers a range from 1 (first script) to the total number of scripts. A value of -1 sets the currently selected script slot.
- **Script file path size** (int) is the number of characters of the script file path string
- **Script file path** (string) is the path of the script file to load
- **Load session** (unsigned int32) automatically loads the scripts from the session file bypassing the script file path argument, where 0=False and 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.Save

Saves the current script in the specified .ini file.

Arguments:

- **Script index** (int) sets the index of the script that will be saved and covers a range from 1 (first script) to the total number of scripts. A value of -1 uses the currently selected script.
- **Script file path size** (int) is the number of characters of the script file path string
- **Script file path** (string) is the path of the script file to save
- **Save session** (unsigned int32) automatically saves the current script into the session file bypassing the script file path argument, where 0=False and 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.Deploy

Deploys a script in the script module.

Arguments:

- **Script index** (int) sets the script to deploy and covers a range from 1 (first script) to the total number of scripts. A value of -1 sets the currently selected script to deploy.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.Undeploy

Undeploys a script in the script module.

Arguments:

- **Script index** (int) sets the script to undeploy and covers a range from 1 (first script) to the total number of scripts. A value of -1 sets the currently selected script to undeploy.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.Run

Runs a script in the script module.

Arguments:

- **Script index** (int) sets the script to run and covers a range from 1 (first script) to the total number of scripts. A value of -1 sets the currently selected script to run.
- **Wait until script finishes** (unsigned int32), where 0=False and 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.Stop

Stops the running script in the script module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.ChsGet

Returns the list of acquired channels in the Script module.

Arguments:

- **Acquire buffer** (unsigned int16) sets the Acquire Buffer number from which to read the list of channels. Valid values are 1 (=Acquire Buffer 1) and 2 (=Acquire Buffer 2).

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of recorded channels. The index is comprised between 0 and 127, and it corresponds to the full list of signals available in the system.
To get the signal name and its corresponding index in the list of the 128 available signals in the Nanonis Controller, use the *Signal.NamesGet* function, or check the RT Idx value in the Signals Manager module.
- **Error** described in the Response message>Body section

Script.ChsSet

Sets the list of acquired channels in the Script module.

Arguments:

- **Acquire buffer** (unsigned int16) sets the Acquire Buffer number from which to set the list of channels. Valid values are 1 (=Acquire Buffer 1) and 2 (=Acquire Buffer 2).
- **Number of channels** (int) is the number of recorded channels. It defines the size of the Channel indexes array
- **Channel indexes** (1D array int) are the indexes of recorded channels. The index is comprised between 0 and 127, and it corresponds to the full list of signals available in the system.
To get the signal name and its corresponding index in the list of the 128 available signals in the Nanonis Controller, use the *Signal.NamesGet* function, or check the RT Idx value in the Signals Manager module.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.DataGet

Returns the data acquired in the Script module.

Arguments:

- **Acquire buffer** (unsigned int16) sets the Acquire Buffer number from which to read the acquired data. Valid values are 1 (=Acquire Buffer 1) and 2 (=Acquire Buffer 2).
- **Sweep number** (int) selects the sweep this function will return the data from. Each sweep is configured as such in the script and it corresponds to each plot displayed in the graphs of the Script module. The sweep numbers start at 0.

Return arguments (if Send response back flag is set to True when sending request message):

- **Data rows** (int) defines the number of rows of the Data array
- **Data columns** (int) defines the number of columns of the Data array
- **Data** (2D array float32) returns the script data
- **Error** described in the Response message>Body section

Script.Autosave

Saves automatically to file the data stored in the Acquire Buffers after running a script in the Script module.

Arguments:

- **Acquire buffer** (unsigned int16) sets the Acquire Buffer number from which to save the data. Valid values are 0 (=Acquire Buffer 1 & Acquire Buffer 2), 1 (=Acquire Buffer 1), and 2 (=Acquire Buffer 2).
- **Sweep number** (int) selects the sweep this function will save the data for. Each sweep is configured as such in the script and it corresponds to each plot displayed in the graphs of the Script module. The sweep numbers start at 0. A value of -1 saves all acquired sweeps.
- **All sweeps to same file** (unsigned int32) decides if all sweeps defined by the Sweep number parameter are saved to the same file (=1) or not (=0).
- **Folder path path size** (int) is the number of characters of the folder path string
- **Folder path** (string) is the folder where the file will be saved. If nothing is sent, the saving routine uses the last used path.
- **Basename path size** (int) is the number of characters of the basename string
- **Basename** (string) is the basename of the file to save. If nothing is sent, the saving routine uses the last used basename.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.LUTOpen

Opens the LUT (Look Up Table) Editor from the Script module.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.LUTLoad

Loads a LUT from file or directly by setting an array of values into the LUT Editor in the script module.

Arguments:

- **LUT index** (int) sets the LUT to load and covers a range from 1 (first LUT) to the total number of LUTs
- **File path size** (int) is the number of characters of the File path string
- **File path** (string) is the path of the file containing a list of values to load into the selected LUT. The extension of the file must be .luts. If nothing is sent, it will use the values sent by the LUT values argument.
- **LUT values size** (int) is the size of the LUT values array
- **LUT values** (1D array float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.LUTSave

Saves a LUT to file.

Arguments:

- **LUT index** (int) sets the LUT to save and covers a range from 1 (first LUT) to the total number of LUTs
- **File path size** (int) is the number of characters of the File path string
- **File path** (string) is the path of the file where the selected LUT will be saved. The extension of the file must be .luts.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Script.LUTDeploy

Deploys a LUT from the LUT Editor in the script module.

Arguments:

- **LUT index** (int) sets the LUT to deploy and covers a range from 1 (first LUT) to the total number of LUTs
- **Wait until finished** (unsigned int32) decides if the function waits until the LUT is fully deployed, as it might take some time.
- **Timeout (ms)** (int) sets the number of milliseconds to wait until the LUT is fully deployed. -1 waits forever.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Lock-In

LockIn.ModOnOffSet

Turns the specified Lock-In modulator on or off.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)
- **Lock-In On/Off** (unsigned int32) turns the specified modulator on or off, where 0=Off and 1=On

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.ModOnOffGet

Returns if the specified Lock-In modulator is turned on or off.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Lock-In On/Off** (unsigned int32) returns if the specified modulator is turned on or off, where 0=Off and 1=On
- **Error** described in the Response message>Body section

LockIn.ModSignalSet

Selects the modulated signal of the specified Lock-In modulator.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)
- **Modulator Signal Index** (int) is the signal index out of the list of 128 signals available in the software. To get a list of the available signals, use the *Signals.NamesGet* function.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.ModSignalGet

Returns the modulated signal of the specified Lock-In modulator.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Modulator Signal Index** (int) is the signal index out of the list of 128 signals available in the software. To get a list of the available signals, use the *Signals.NamesGet* function
- **Error** described in the Response message>Body section

LockIn.ModPhasRegSet

Sets the phase register index of the specified Lock-In modulator.

Each modulator can work on any phase register (frequency). Use this function to assign the modulator to one of the 8 available phase registers (index 1-8).

Use the *LockIn.ModPhaFreqSet* function to set the frequency of the phase registers.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)
- **Phase Register Index** (int) is the index of the phase register of the specified Lock-In modulator. Valid values are index 1 to 8.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.ModPhasRegGet

Returns the phase register index of the specified Lock-In modulator.

Each modulator can work on any phase register (frequency generator).

Use the *LockIn.ModPhaseRegFreqGet* function to get the frequency of the phase registers.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Phase Register Index** (int) is the index of the phase register of the specified Lock-In modulator. Valid values are index 1 to 8
- **Error** described in the Response message>Body section

LockIn.ModHarmonicSet

Sets the harmonic of the specified Lock-In modulator.

The modulator is bound to a phase register (frequency generator), but it can work on harmonics. Harmonic 1 is the base frequency (the frequency of the frequency generator).

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)
- **Harmonic** (int) is the harmonic of the specified Lock-In modulator. Valid values start from 1 (=base frequency)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.ModHarmonicGet

Returns the harmonic of the specified Lock-In modulator.

The modulator is bound to a phase register (frequency generator), but it can work on harmonics. Harmonic 1 is the base frequency (the frequency of the frequency generator).

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Harmonic** (int) is the harmonic of the specified Lock-In modulator. Valid values start from 1 (=base frequency)
- **Error** described in the Response message>Body section

LockIn.ModPhasSet

Sets the modulation phase offset of the specified Lock-In modulator.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)
- **Phase (deg)** (float32) is the modulation phase offset of the specified Lock-In modulator

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.ModPhasGet

Returns the modulation phase offset of the specified Lock-In modulator.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Phase (deg)** (float32) is the modulation phase offset of the specified Lock-In modulator
- **Error** described in the Response message>Body section

LockIn.ModAmpSet

Sets the modulation amplitude of the specified Lock-In modulator.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)
- **Amplitude** (float32) is the modulation amplitude of the specified Lock-In modulator

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.ModAmpGet

Returns the modulation amplitude of the specified Lock-In modulator.

Arguments:

- **Modulator number** (int) is the number that specifies which modulator to use. It starts from number 1 (=Modulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Amplitude** (float32) is the modulation amplitude of the specified Lock-In modulator
- **Error** described in the Response message>Body section

LockIn.ModPhasFreqSet

Sets the frequency of the specified Lock-In phase register/modulator.

The Lock-in module has a total of 8 frequency generators / phase registers. Each modulator and demodulator can be bound to one of the phase registers.

This function sets the frequency of one of the phase registers.

Arguments:

- **Modulator number** (int) is the number that specifies which phase register/modulator to use. It starts from number 1 (=Modulator 1)
- **Frequency (Hz)** (float64) is the frequency of the specified Lock-In phase register

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.ModPhasFreqGet

Returns the frequency of the specified Lock-In phase register/modulator.

The Lock-in module has a total of 8 frequency generators / phase registers. Each modulator and demodulator can be bound to one of the phase registers.

This function gets the frequency of one of the phase registers.

Arguments:

- **Modulator number** (int) is the number that specifies which phase register/modulator to use. It starts from number 1 (=Modulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Frequency (Hz)** (float64) is the frequency of the specified Lock-In phase register
- **Error** described in the Response message>Body section

LockIn.DemodSignalSet

Selects the demodulated signal of the specified Lock-In demodulator.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)
- **Demodulator Signal Index** (int) is the signal index out of the list of 128 signals available in the software. To get a list of the available signals, use the *Signals.NamesGet* function.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.DemodSignalGet

Returns the demodulated signal of the specified Lock-In demodulator.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Demodulator Signal Index** (int) is the signal index out of the list of 128 signals available in the software. To get a list of the available signals, use the *Signals.NamesGet* function
- **Error** described in the Response message>Body section

LockIn.DemodHarmonicSet

Sets the harmonic of the specified Lock-In demodulator.

The demodulator demodulates the input signal at the specified harmonic overtone of the frequency generator. Harmonic 1 is the base frequency (the frequency of the frequency generator).

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)
- **Harmonic** (int) is the harmonic of the specified Lock-In demodulator. Valid values start from 1 (=base frequency)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.DemodHarmonicGet

Returns the harmonic of the specified Lock-In demodulator.

The demodulator demodulates the input signal at the specified harmonic overtone of the frequency generator. Harmonic 1 is the base frequency (the frequency of the frequency generator).

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Harmonic** (int) is the harmonic of the specified Lock-In demodulator. Valid values start from 1 (=base frequency)
- **Error** described in the Response message>Body section

LockIn.DemodHPFilterSet

Sets the properties of the high-pass filter applied to the demodulated signal of the specified demodulator.

The high-pass filter is applied on the demodulated signal before the actual demodulation. It is used to get rid of DC or low-frequency components which could result in undesired components close to the modulation frequency on the demodulator output signals (X,Y).

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)
- **HP Filter Order** (int) is the high-pass filter order. Valid values are from -1 to 8, where -1=no change, 0=filter off.
- **HP Filter Cutoff Frequency (Hz)** (float32) is the high-pass filter cutoff frequency in Hz, where 0 = no change.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.DemodHPFilterGet

Returns the properties of the high-pass filter applied to the demodulated signal of the specified demodulator.

The high-pass filter is applied on the demodulated signal before the actual demodulation. It is used to get rid of DC or low-frequency components which could result in undesired components close to the modulation frequency on the demodulator output signals (X,Y).

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **HP Filter Order** (int) is the high-pass filter order. Valid values are from 0 to 8, where 0=filter off
- **HP Filter Cutoff Frequency (Hz)** (float32) is the high-pass filter cutoff frequency in Hz
- **Error** described in the Response message>Body section

LockIn.DemodLPFilterSet

Sets the properties of the low-pass filter applied to the demodulated signal of the specified demodulator.

The low-pass filter is applied on the demodulator output signals (X,Y) to remove undesired components. Lower cut-off frequency means better suppression of undesired frequency components, but longer response time (time constant) of the filter.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)
- **LP Filter Order** (int) is the low-pass filter order. Valid values are from -1 to 8, where -1=no change, 0=filter off.
- **LP Filter Cutoff Frequency (Hz)** (float32) is the low-pass filter cutoff frequency in Hz, where 0 = no change.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.DemodLPFilterGet

Returns the properties of the low-pass filter applied to the demodulated signal of the specified demodulator.

The low-pass filter is applied on the demodulator output signals (X,Y) to remove undesired components. Lower cut-off frequency means better suppression of undesired frequency components, but longer response time (time constant) of the filter.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **LP Filter Order** (int) is the low-pass filter order. Valid values are from -1 to 8, where -1=no change, 0=filter off.
- **LP Filter Cutoff Frequency (Hz)** (float32) is the low-pass filter cutoff frequency in Hz, where 0 = no change.
- **Error** described in the Response message>Body section

LockIn.DemodPhasRegSet

Sets the phase register index of the specified Lock-In demodulator.

Each demodulator can work on any phase register (frequency). Use this function to assign the demodulator to one of the 8 available phase registers (index 1-8).

Use the *LockIn.ModPhaFreqSet* function to set the frequency of the phase registers.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)
- **Phase Register Index** (int) is the index of the phase register of the specified Lock-In demodulator. Valid values are index 1 to 8.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.DemodPhasRegGet

Returns the phase register index of the specified Lock-In demodulator.

Each demodulator can work on any phase register (frequency). Use the *LockIn.ModPhaFreqSet* function to set the frequency of the phase registers.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Phase Register Index** (int) is the index of the phase register of the specified Lock-In demodulator. Valid values are index 1 to 8.
- **Error** described in the Response message>Body section

LockIn.DemodPhasSet

Sets the reference phase of the specified Lock-In demodulator.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)
- **Phase (deg)** (float32) is the reference phase of the specified Lock-In demodulator

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.DemodPhasGet

Returns the reference phase of the specified Lock-In demodulator.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Phase (deg)** (float32) is the reference phase of the specified Lock-In demodulator
- **Error** described in the Response message>Body section

LockIn.DemodSyncFilterSet

Switches the synchronous (Sync) filter of the specified demodulator On or Off.

The synchronous filter is applied on the demodulator output signals (X,Y) after the low-pass filter. It is very good in suppressing harmonic components (harmonics of the demodulation frequency), but does not suppress other frequencies.

The sync filter does not output a continuous signal, it only updates the value after each period of the demodulation frequency.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)
- **Sync Filter** (unsigned int32) switches the synchronous filter of the specified demodulator on or off, where 0=Off and 1=On

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.DemodSyncFilterGet

Returns the status (on/off) of the synchronous (Sync) filter of the specified demodulator.

The synchronous filter is applied on the demodulator output signals (X,Y) after the low-pass filter. It is very good in suppressing harmonic components (harmonics of the demodulation frequency), but does not suppress other frequencies.

The sync filter does not output a continuous signal, it only updates the value after each period of the demodulation frequency.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Sync Filter** (unsigned int32) is the synchronous filter of the specified demodulator, where 0=Off and 1=On
- **Error** described in the Response message>Body section

LockIn.DemodRTSignalsSet

Sets the signals available for acquisition on the real-time system from the specified demodulator.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)
- **RT Signals** (unsigned int32) sets which signals from the specified demodulator should be available on the Real-time system. 0 sets the available RT Signals to X/Y, 1 sets them to R/phi.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockIn.DemodRTSignalsGet

Returns which the signals are available for acquisition on the real-time system from the specified demodulator.

Arguments:

- **Demodulator number** (int) is the number that specifies which demodulator to use. It starts from number 1 (=Demodulator 1)

Return arguments (if Send response back flag is set to True when sending request message):

- **RT Signals** (unsigned int32) returns which signals from the specified demodulator are available on the Real-time system. 0 means X/Y, and 1 means R/phi.
- **Error** described in the Response message>Body section

Lock-In Frequency Sweep

LockInFreqSwp.Open

Opens the Transfer function (Lock-In Frequency Sweep) module.

The transfer function does not run when its front panel is closed. To automate measurements it might be required to open the module first using this VI.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockInFreqSwp.Start

Starts a Lock-In frequency sweep.

Arguments:

- **Get Data** (unsigned int32) defines if the function returns the recorder channels and data
- **Direction** (unsigned int32) sets the direction of the frequency sweep. 0 means sweep down (from upper limit to lower limit) and 1 means sweep up (from lower limit to upper limit)

Return arguments (if Send response back flag is set to True when sending request message):

- **Channels names size** (int) is the size in bytes of the recorder channels names array
- **Channels names number** (int) is the number of elements of the recorded channels names array
- **Channels names** (1D array string) returns the array of recorded channel names (strings), where each string comes prepended by its size in bytes
- **Data rows** (int) is the number of rows of the returned data array (the first row is the swept frequency, and each additional row contains the data of each recorded channel)
- **Data columns** (int) is the number of recorded points (number of steps plus 1)
- **Data** (2D array float32) returns the recorded data. The number of rows is defined by *Data rows*, and the number of columns is defined by *Data columns*
- **Error** described in the Response message>Body section

LockInFreqSwp.SignalSet

Sets the sweep signal used in the Lock-In frequency sweep module.

Arguments:

- **Sweep signal index** (int) sets the sweep signal index out of the list of sweep signals to use, where -1 means no signal selected

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockInFreqSwp.SignalGet

Returns the sweep signal used in the Lock-In frequency sweep module.

Arguments:

Return arguments (if Send response back flag is set to True when sending request message):

- **Sweep signal index** (int) is the sweep signal index selected out of the list of sweep signals, where -1 means no signal selected
- **Error** described in the Response message>Body section

LockInFreqSwp.LimitsSet

Sets the frequency limits in the Lock-In frequency sweep module.

Arguments:

- **Lower limit (Hz)** (float32) sets the lower frequency limit in Hz
- **Upper limit (Hz)** (float32) sets the lower frequency limit in Hz

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockInFreqSwp.LimitsGet

Returns the frequency limits in the Lock-In frequency sweep module.

Arguments:

Return arguments (if Send response back flag is set to True when sending request message):

- **Lower limit (Hz)** (float32) is the lower frequency limit in Hz
- **Upper limit (Hz)** (float32) is the lower frequency limit in Hz
- **Error** described in the Response message>Body section

LockInFreqSwp.PropsSet

Sets the configuration of the Transfer Function (Lock-In frequency sweep) module.

Arguments:

- **Number of steps** (unsigned int16) is the number of frequency steps over the sweep range (logarithmic distribution). The number of data points = number of steps + 1. If set to 0, the number of steps is left unchanged
- **Integration periods** (unsigned int16) is the number of Lock in periods to average for one measurement.
- **Minimum integration time (s)** (float32) is the minimum integration time in seconds to average each measurement
- **Settling periods** (unsigned int16) is the number of Lock in periods to wait before acquiring data at each point of the sweep
- **Minimum Settling time (s)** (float32) is the minimum settling time in seconds to wait before acquiring data at each point of the sweep
- **Autosave** (unsigned int32) automatically saves the data at end of sweep
- **Save dialog** (unsigned int32) will open a dialog box when saving the data
- **Basename size** (int) is the size (number of characters) of the basename string
- **Basename** (string) is the basename of the saved files

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

LockInFreqSwp.PropsGet

Returns the configuration of the Transfer Function (Lock-In frequency sweep) module.

Arguments:

Return arguments (if Send response back flag is set to True when sending request message):

- **Number of steps** (unsigned int16) is the number of frequency steps over the sweep range (logarithmic distribution). The number of data points = number of steps + 1
- **Integration periods** (unsigned int16) is the number of Lock in periods to average for one measurement.
- **Minimum integration time (s)** (float32) is the minimum integration time in seconds to average each measurement
- **Settling periods** (unsigned int16) is the number of Lock in periods to wait before acquiring data at each point of the sweep
- **Minimum Settling time (s)** (float32) is the minimum settling time in seconds to wait before acquiring data at each point of the sweep
- **Autosave** (unsigned int32) automatically saves the data at end of sweep
- **Save dialog** (unsigned int32) will open a dialog box when saving the data
- **Basename size** (int) is the size (number of characters) of the basename string
- **Basename** (string) is the basename of the saved files
- **Error** described in the Response message>Body section

Utilities

Util.SessionPathGet

Returns the session path.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Session path size** (int) is the number of characters of the Session path string
- **Session path** (string)
- **Error** described in the Response message>Body section

Util.SessionPathSet

Sets the session folder path.

Arguments:

- **Session path size** (int) is the number of characters of the Session path string
- **Session path** (string)
- **Save settings to previous** (unsigned int32) determines if the settings are saved to the previous session file before changing it, where 0=False and 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.SettingsLoad

Loads the settings from the specified .ini file.

Arguments:

- **Settings file path size** (int) is the number of characters of the Settings file path string
- **Settings file path** (string) is the path of the settings file to load
- **Load session settings** (unsigned int32) automatically loads the current settings from the session file bypassing the settings file path argument, where 0=False and 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.SettingsSave

Saves the current settings in the specified .ini file.

Arguments:

- **Settings file path size** (int) is the number of characters of the Settings file path string
- **Settings file path** (string) is the path of the settings file to save
- **Save session settings** (unsigned int32) automatically saves the current settings into the session file bypassing the settings file path argument, where 0=False and 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.LayoutLoad

Loads a layout from the specified .ini file.

Arguments:

- **Layout file path size** (int) is the number of characters of the layout file path string
- **Layout file path** (string) is the path of the layout file to load
- **Load session layout** (unsigned int32) automatically loads the layout from the session file bypassing the layout file path argument, where 0=False and 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.LayoutSave

Saves the current layout in the specified .ini file.

Arguments:

- **Layout file path size** (int) is the number of characters of the layout file path string
- **Layout file path** (string) is the path of the layout file to save
- **Save session layout** (unsigned int32) automatically saves the current layout into the session file bypassing the layout file path argument, where 0=False and 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.Lock

Locks the Nanonis software.

Launches the Lock modal window, preventing the user to interact with the Nanonis software until unlocking it manually or through the *Util.Unlock* function.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.Unlock

Unlocks the Nanonis software.

Closes the Lock modal window which prevents the user to interact with the Nanonis software.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

Error described in the Response message>Body section

Util.RTFreqSet

Sets the Real Time controller frequency.

Arguments:

- **RT frequency** (float32) is the Real Time frequency in Hz

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.RTFreqGet

Gets the Real Time controller frequency.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **RT frequency** (float32) is the Real Time frequency in Hz
- **Error** described in the Response message>Body section

Util.AcqPeriodSet

Sets the Acquisition Period (s) in the TCP Receiver.

Arguments:

- **Acquisition Period (s)** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.AcqPeriodGet

Gets the Acquisition Period (s) in the TCP Receiver.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Acquisition Period (s)** (float32)
- **Error** described in the Response message>Body section

Util.RTOversamplSet

Sets the Real-time oversampling in the TCP Receiver.

The 24 signals are oversampled on the RT engine before they are sent to the host. The oversampling affects the maximum Spectrum Analyzer frequency and other displays.

Arguments:

- **RT oversampling** (int)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.RTOversamplGet

Returns the Real-time oversampling in the TCP Receiver.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **RT oversampling** (int)
- **Error** described in the Response message>Body section

Util.Quit

Quits the Nanonis software with the option to save settings, layout and signals (same functionality provided by the dialog window that pops-up when quitting the software through the File menu).

Arguments:

- **Use the stored values** (unsigned int32) automatically ignores the rest of the arguments (0=False and 1=True) and saves settings, layout, and signals according to the last time the software quit. This configuration is stored in the Main-Options settings.ini file located in the Certificate folder.
- **Settings name size** (int) is the number of characters of the Settings name string
- **Settings name** (string) is the name of the settings file to save when quitting. The list of settings can be found in the Settings section of the Main Options under the File menu. If left empty, no settings are saved (unless the argument “Use the stored values” is 1).
- **Layout name size** (int) is the number of characters of the Layout name string
- **Layout name** (string) is the name of the layout file to save when quitting. The list of layouts can be found in the Layouts section of the Main Options under the File menu. If left empty, no layout is saved (unless the argument “Use the stored values” is 1).
- **Save signals** (unsigned int32) automatically saves (0=False and 1=True) the signal configuration currently set in the Signals Manager if it has been changed. The signals configuration is stored in the Signals settings.ini file located in the Certificate folder.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

Util.VersionGet

Returns the version information of the Nanonis software.

Arguments: None

Return arguments (if Send response back flag is set to True when sending request message):

- **Product Line size** (int) is the number of characters of the Product Line string
- **Product Line** (string) returns “Nanonis SPM Control Software” or “Nanonis Tramea Software”.
- **Version size** (int) is the number of characters of the Version string
- **Version** (string) returns the software version (e.g. Generic 5)
- **Host App. Release** (unsigned int32) returns the host application release number.
- **RT Engine release** (unsigned int32) returns the RT Engine application release number.
- **Error** described in the Response message>Body section

File

File.datLoad

Returns the contents (Channel names, Data, and Header) of a .dat file saved with the Nanonis software.

Arguments:

- **File path size** (int) is the number of characters of the File path string
- **File path** (string) is the path of the .dat file to load
- **Read only header** (unsigned int32) defines if only the Header is returned or if also the Channel names & Data are returned, where 0=False and 1=True

Return arguments (if Send response back flag is set to True when sending request message):

- **Channels names size** (int) is the size in bytes of the Channels names string array
- **Number of channels** (int) is the number of elements of the Channels names string array
- **Channels names** (1D array string) returns the list of channels names. The size of each string item comes right before it as integer 32
- **Data rows** (int) defines the number of rows of the Data array
- **Data columns** (int) defines the number of columns of the Data array
- **Data** (2D array float32) returns the data stored in the file
- **Header rows** (int) defines the number of rows of the Header array
- **Header columns** (int) defines the number of columns of the Header array
- **Header** (2D array string) returns the header line by line and the attribute followed by its value (i.e. Attribute 1->Value->Attribute 2->Value...). The size of each string item comes right before it as integer 32
- **Error** described in the Response message>Body section

MCVA5 Preamplifier

MCVA5.UserInSet

Assigns a user input to one of the four channels of the selected preamplifier.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Channel Nr.** (unsigned int16) is the channel number out of the four available channels in the MCVA5 preamplifier, where valid values are 1, 2, 3, 4
- **User Input** (unsigned int32) is the user input number that will be assigned to the selected channel of the selected preamplifier

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

MCVA5.UserInGet

Returns the user input assigned to the selected channel of the selected preamplifier.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Channel Nr.** (unsigned int16) is the channel number out of the four available channels in the MCVA5 preamplifier, where valid values are 1, 2, 3, 4

Return arguments (if Send response back flag is set to True when sending request message):

- **User Input** (unsigned int32) is the user input number that is currently assigned to the selected channel of the selected preamplifier
- **Error** described in the Response message>Body section

MCVA5.GainSet

Sets the gain of one of the four channels of the selected preamplifier.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Channel Nr.** (unsigned int16) is the channel number out of the four available channels in the MCVA5 preamplifier, where valid values are 1, 2, 3, 4
- **Gain** (unsigned int16) is the gain that will be applied to the selected channel of the selected preamplifier, where 0=1, 1=10, 2=100, 3=1000

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

MCVA5.GainGet

Returns the gain applied to one of the four channels of the selected preamplifier.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Channel Nr.** (unsigned int16) is the channel number out of the four available channels in the MCVA5 preamplifier, where valid values are 1, 2, 3, 4

Return arguments (if Send response back flag is set to True when sending request message):

- **Gain** (unsigned int16) is the gain currently applied to the selected channel of the selected preamplifier, where 0=1, 1=10, 2=100, 3=1000
- **Error** described in the Response message>Body section

MCVA5.InputModeSet

Sets the input mode of one of the four channels of the selected preamplifier.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Channel Nr.** (unsigned int16) is the channel number out of the four available channels in the MCVA5 preamplifier, where valid values are 1, 2, 3, 4
- **Input mode** (unsigned int16) is the input mode that will be applied to the selected channel of the selected preamplifier, where 0=A-B, 1=A, 2=GND

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

MCVA5.InputModeGet

Returns the input mode applied to one of the four channels of the selected preamplifier.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Channel Nr.** (unsigned int16) is the channel number out of the four available channels in the MCVA5 preamplifier, where valid values are 1, 2, 3, 4

Return arguments (if Send response back flag is set to True when sending request message):

- **Input mode** (unsigned int16) is the input mode currently applied to the selected channel of the selected preamplifier, where 0=A-B, 1=A, 2=GND
- **Error** described in the Response message>Body section

MCVA5.CouplingSet

Sets the coupling mode of one of the four channels of the selected preamplifier.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Channel Nr.** (unsigned int16) is the channel number out of the four available channels in the MCVA5 preamplifier, where valid values are 1, 2, 3, 4
- **Coupling** (unsigned int16) is the coupling mode that will be applied to the selected channel of the selected preamplifier, where 0=DC and 1=AC

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

MCVA5.CouplingGet

Returns the coupling mode of one of the four channels of the selected preamplifier.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Channel Nr.** (unsigned int16) is the channel number out of the four available channels in the MCVA5 preamplifier, where valid values are 1, 2, 3, 4

Return arguments (if Send response back flag is set to True when sending request message):

- **Coupling** (unsigned int16) is the coupling mode that will be applied to the selected channel of the selected preamplifier, where 0=DC and 1=AC
- **Error** described in the Response message>Body section

MCVA5.SingleStateUpdate

Reads once the state of the selected preamplifier and returns the Ready and Overload states of all channels.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4

Return arguments (if Send response back flag is set to True when sending request message):

- **Channel 1 Ready** (unsigned int32) where 0=no ready, 1=ready
- **Channel 1 Overload** (unsigned int32) where 0=not overload, 1=overload
- **Channel 2 Ready** (unsigned int32) where 0=no ready, 1=ready
- **Channel 2 Overload** (unsigned int32) where 0=not overload, 1=overload
- **Channel 3 Ready** (unsigned int32) where 0=no ready, 1=ready
- **Channel 3 Overload** (unsigned int32) where 0=not overload, 1=overload
- **Channel 4 Ready** (unsigned int32) where 0=no ready, 1=ready
- **Channel 4 Overload** (unsigned int32) where 0=not overload, 1=overload
- **Error** described in the Response message>Body section

MCVA5.ContStateUpdateSet

Configures the selected preamplifier to read continuously the state of all channels.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Continuous Read** (unsigned int32) where 0=false, 1=true

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

MCVA5.ContStateUpdateGet

Returns if the selected preamplifier is configured to read continuously the state of all channels.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4

Return arguments (if Send response back flag is set to True when sending request message):

- **Continuous Read** (unsigned int32) where 0=false, 1=true
- **Error** described in the Response message>Body section

MCVA5.SingleTempUpdate

Reads once the temperatures of all channels of the selected preamplifier.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4

Return arguments (if Send response back flag is set to True when sending request message):

- **Channel 1** (unsigned int16) is the temperature in °C of Channel 1
- **Channel 2** (unsigned int16) is the temperature in °C of Channel 2
- **Channel 3** (unsigned int16) is the temperature in °C of Channel 3
- **Channel 4** (unsigned int16) is the temperature in °C of Channel 4
- **Error** described in the Response message>Body section

MCVA5.ContTempUpdateSet

Configures the selected preamplifier to read continuously the temperature of all channels.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4
- **Continuous Read** (unsigned int32) where 0=false, 1=true

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

MCVA5.ContTempUpdateGet

Returns if the selected preamplifier is configured to read continuously the temperature of all channels.

Arguments:

- **Preamplifier Nr.** (unsigned int16) is the preamplifier module number, where valid values are 1, 2, 3, 4

Return arguments (if Send response back flag is set to True when sending request message):

- **Continuous Read** (unsigned int32) where 0=false, 1=true
- **Error** described in the Response message>Body section

Generic PI Controller

PICtrl.OnOffSet

Switches the Generic PI Controller On or Off.

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.
- **Controller status** (unsigned int32) switches the controller Off (=0) or On (=1)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

PICtrl.OnOffGet

Returns the status of the Generic PI Controller..

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.

Return arguments (if Send response back flag is set to True when sending request message):

- **Controller status** (unsigned int32) indicates if the controller is Off (=0) or On (=1)
- **Error** described in the Response message>Body section

PICtrl.CtrlChSet

Sets the index of the output controlled by the Generic PI controller.

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.
- **Control signal index** (int) sets the output index to be used.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

PICtrl.CtrlChGet

Gets the index of the output controlled by the Generic PI controller.

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.

Return arguments (if Send response back flag is set to True when sending request message):

- **Control signal index** (int) returns the output index to be used.
- **Control signals names size** (int) is the size in bytes of the Control Signals Names array
- **Number of Control signals names** (int) is the number of elements of the Control Signals Names array
- **Control signals names** (1D array string) returns an array of Control Signals Names. Each element of the array is preceded by its size in bytes
- **Number of Control signals indexes** (int) is the number of elements of the Control Signals Indexes array
- **Control signals indexes** (1D array int) returns an array of Control Signals Indexes
- **Error** described in the Response message>Body section

PICtrl.InputChSet

Sets the index of the input channel in the Generic PI controller.

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.
- **Input index** (int) sets the input index to be used.

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

PICtrl.InputChGet

Gets the index of the input channel in the Generic PI controller.

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.

Return arguments (if Send response back flag is set to True when sending request message):

- **Input index** (int) returns the input index to be used.
- **Input signals names size** (int) is the size in bytes of the Input Signals Names array
- **Number of Input signals names** (int) is the number of elements of the Input Signals Names array
- **Input signals names** (1D array string) returns an array of Input Signals Names. Each element of the array is preceded by its size in bytes
- **Number of Input signals indexes** (int) is the number of elements of the Input Signals Indexes array
- **Input signals indexes** (1D array int) returns an array of Input Signals Indexes
- **Error** described in the Response message>Body section

PICtrl.PropsSet

Sets the properties of the Generic PI controller.

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.
- **Setpoint** (float32)
- **P gain** (float32)
- **I gain** (float32)
- **Slope** (unsigned int16) where 0 means no change, 1 means Positive, and 2 means Negative

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

PICtrl.PropsGet

Gets the properties of the Generic PI controller.

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.

Return arguments (if Send response back flag is set to True when sending request message):

- **Setpoint** (float32)
- **P gain** (float32)
- **I gain** (float32)
- **Slope** (unsigned int16) where 0 means Positive, and 1 means Negative
- **Error** described in the Response message>Body section

PICtrl.CtrlChPropsSet

Sets the properties of the control signal in the Generic PI controller.

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.
- **Lower limit** (float32)
- **Upper limit** (float32)

Return arguments (if Send response back flag is set to True when sending request message):

- **Error** described in the Response message>Body section

PICtrl.CtrlChPropsGet

Gets the properties of the control signal in the Generic PI controller.

Arguments:

- **Controller index** (int) is the index of the controller. Valid values are 1 to 8.

Return arguments (if Send response back flag is set to True when sending request message):

- **Lower limit** (float32)
- **Upper limit** (float32)
- **Error** described in the Response message>Body section